

C

## COMPLETE CATALOG



System provider for Incontinence Care and more.





*With our range of draining incontinence products, we at GHC German Health Care GmbH have focused on a care area in which we can offer a comprehensive product range with a very high level of expertise. And our customers can rely on this, especially during challenging times.*

## Hello from Flensburg

You will find a wide selection of transurethral and suprapubic balloon catheters, various catheter insertion kits, male external catheters (MEC), a wide range of urine collecting systems and related accessories. We are also excited to introduce the first versions of our new IC product line, CARE CATH Liv & Leiv (pages 62/63), along with intermittent catheters from Curan (pp. 64–65).

In our current catalogue, you will discover additional new products that we have added for our customers in home-care providers, medical supply retailers, and hospitals. Our CARE CON condom-catheter series (pp. 38–40) offers excellent, reliable quality at an attractive price point.

New to our range are the MicroN Safety safety port needles (pp. 66–67) for parenteral therapy, pain management, and oncology infusion treatment, as well as GHC's paediatric line-up—including the 1 000 ml drip chamber pouch, the 210 ml leg bag, and a Velcro fastener set (pp. 52–53). In addition to our stringent product-quality standards, we guarantee reliable delivery for these new items.

Our comprehensive catalogue gives you an overview of the full GHC range. For the latest information and newly launched products, please visit [www.ghc-medical.de](http://www.ghc-medical.de).

Thanks to expanded warehouse capacity, proactive procurement planning, and a strong service focus—complete with dedicated support offerings—we are well equipped for the future. For example, we have developed a customer portal in which all MDR-relevant documentation for our product range is accessible and available for download at any time. This service saves considerable time and ensures you always have current certificates and technical files on hand.

Please feel free to contact us to discover how a partnership with GHC German Health Care GmbH can benefit your organization.

We look forward to speaking with you.

Nordic greetings from the GHC team

# C

Introduction	page 3
Table of Contents	page 4/5
Quality policy of GHC	page 6
Catheter guide	page 7

# C

## Draining incontinence

### CARE FLOW transurethral catheters

CARE FLOW silicone catheter, Nelaton tip	page 8
CARE FLOW Soft silicone catheter, Nelaton tip	page 9
CARE FLOW silicone catheter - <i>integral</i> , Nelaton tip, with spigot	page 10
BLUECATH silicone catheter - <i>integral</i> , Nelaton tip, with spigot	page 11
CARE FLOW silicone elastomer catheter, Nelaton tip	page 12
CARE FLOW silicone catheter, Tiemann tip	page 13
CARE FLOW silicone catheter, Flute tip	page 14
CARE FLOW silicone elastomer catheter, Couvelaire tip	page 15
CARE FLOW Siliconised latex catheters, Nelaton tip	page 16
CARE FLOW PVC disposable catheter, Nelaton & Tiemann tip	page 17

### CARE FLOW irrigation catheter

CARE FLOW silicone irrigation catheter Nelaton tip	page 18
CARE FLOW silicone irrigation catheter Dufour tip	page 19
CARE FLOW silicone irrigation catheter Couvelaire tip	page 20
CARE FLOW silicone irrigation catheter Mercier tip	page 21

### CARE FLOW suprapubic catheters

CARE FLOW UNIBAL silicone catheter - <i>integral</i> - centre open	page 22
CARE FLOW UNIBAL silicone catheter - <i>integral</i> - centre open, with spigot	page 23
CARE FLOW silicone catheter - <i>integral</i> - centre open, with spigot	page 24
CARE FLOW UNIBAL silicone catheter - <i>integral</i> - centre open, short	page 25
BLUECATH silicone catheter - <i>integral</i> - dilatation tip, with spigot	page 26
BLUECATH silicone catheter - <i>integral</i> - centre open, with spigot	page 27
CARE FLOW silicone catheter <i>blue - integral</i> - centre open, with spigot	page 28
CARE FLOW silicone catheter - <i>integral</i> - centre open, short, with spigot	page 29
CARE FLOW UNIBAL puncture set	page 30
CARE FLOW UNIBAL replacement set, incl. PTFE guide wire	page 31

### CARE FLOW nephrostomy catheter

CARE FLOW silicone catheter - <i>integral</i> - Dilatation tip	page 32
--	---------

### CARE SET Catheter insertion kits and accessories

Catheter insertion kits Basic and Clinical	page 33
Catheter insertion kit Comfort and Comfort Plus	page 34
Catheter insertion kits Comfort Eco	page 35
Accessories (blocker syringe, lubricant gel, catheter spigot, velcro fixing plaster)	page 36

### CARE FLOW catheter valves

CARE FLOW and CARE VENT	page 37
-------------------------	---------

# C

### CARE CON silicone male external catheters (mec)

CARE CON male external catheters (MEC) comfort /Type I	page 38
CARE CON male external catheters (MEC) special /Type II	page 39
CARE CON male external catheters (MEC) active / type III	page 40

### CARE FLOW urine drainage systems and accessories

CARE FLOW UB 4000TUR Closed irrigation bag system	page 41
CARE FLOW UB 4000 Closed urine drainage system	page 42
CARE FLOW UB 2000 Closed urine drainage system	page 43
CARE FLOW UB 2000S Closed urine drainage system, CARE BAG	page 44
CARE FLOW UB 2000A Bedside bag	page 45
CARE FLOW UB 2000E, UB 2000US Disposable bed bags, PVC bag holder	page 46
CARE FLOW UB 2000USo, UB 1500USo, SB 1500USo Disposable bag without drain	page 47
CARE FLOW UB 1000 Closed urine drainage system day & night	page 48
CARE FLOW UB 750/750A/750US 1-chamber leg bag system	page 49
CARE FLOW UB 750F/750UF 1-chamber leg bag system with flip valve	page 50
CARE FLOW UB 500/500US 3-chamber leg bag system	page 51

### CARE FLOW Paediatric Care

CARE FLOW Closed urine drainage system for children <b>new</b>	page 52
CARE FLOW Leg Bag Velcro Strap Set for Children <b>new</b>	page 53

### CARE FIX leg bag fixation systems

CARE FIX Universal Velcro strap	page 54
CARE FIX Velcro straps premium set, with and without buttons	page 55
CARE FIX Fixation stockings with integrated leg bag pocket	page 56
CARE FIX Leg warmers for upper and lower leg cuffs	page 57
CARE FIX Premium one leg pants	page 58

### CARE DRESS dressings

Suprapubic dressing set	page 59
PEG dressing set	page 60

# C

## Intermittent catheterisation (IC)

### CARE IK Intermittent catheterisation

CARE CATH <sup>IK</sup> IC disposable catheter, hydrophilic, 17 cm	page 62
CARE CATH <sup>IK</sup> IC disposable catheter, hydrophilic, 41 cm	page 63
Curan Lady & Man	page 64
Curan Complete & Closed	page 65

# C

## Oncology

MicroN Safety port cannula	page 66/67
----------------------------	------------

# C



*German Health Care GmbH stands for innovative solutions in the medical care for people with physical disabilities. As an importer and distributor of products specifically for bladder management, we offer our customers a broad portfolio of customised solutions to improve their customers' quality of life - with the assurance of high product quality and reliable use.*

We have developed our extensive product portfolio in strategic partnership with renowned manufacturers. Through ongoing dialogue, our items are continuously refined, and new medical devices are regularly added. Market-leading home-care providers, hospitals, and purchasing cooperatives in the medical-supply sector have relied on the quality of our products and services for years.


All manufacturing is conducted according to precise specifications across our partners' global production sites. We routinely audit these facilities to verify the integrity of their manufacturing processes. After production, every batch undergoes comprehensive quality control by the manufacturer, and once again upon arrival at our high-bay warehouse in Wees. Only flawlessly produced items are released for distribution. This rigorous quality policy is ingrained in our operations—our customers can depend on it.

Moreover, our Quality Management System, certified to DIN EN ISO 13485, is regularly monitored and audited by DQS.

Equally fundamental to our success are our people. Our qualified in-house and field-based medical-device consultants are at your disposal for any questions regarding the GHC range and ensure that your inquiries are handled and dispatched as swiftly as possible. Continuous in-house workshops, ongoing professional development, and IT-supported processes guarantee this expertise, so you can rest assured that your orders and deliveries will be processed without delay.

By partnering with GHC German Health Care GmbH, you secure top-tier quality in incontinent care solutions and beyond.

# C

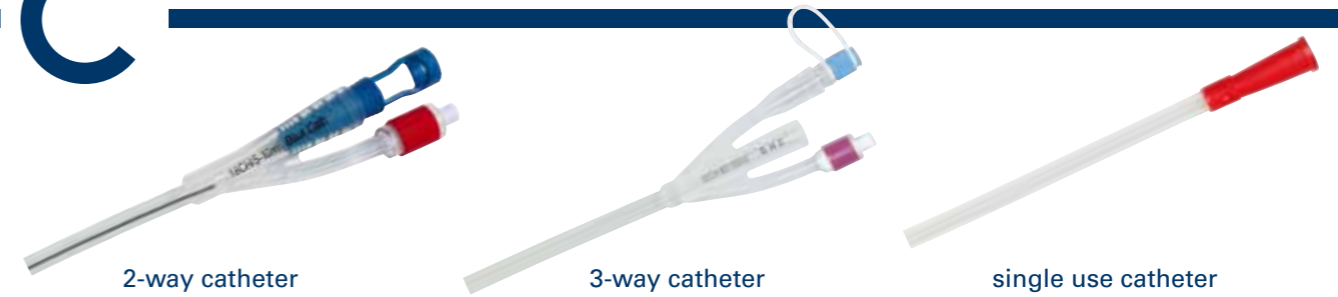

2-way catheter

3-way catheter

single use catheter


1 Charriere (CH) = 1 French (FR) = 0,33mm

We carry catheters of different materials with different specifications to meet all the needs of our customers. Among them are catheters made of full silicone and Silicone elastomer for long-term drainage, but also latex catheter for short-term drainage.

### Classic catheter materials of catheters for long-term drainage



100% Silicone

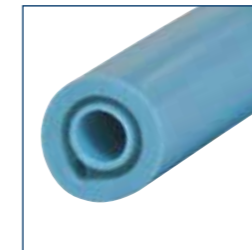

Silicone elastomer – Latex core with silicone coating



only for short-term usage, up to a maximum of 5 days

### Catheter tips for transurethral and suprapubic urine drainage



Nelaton tip



Tiemann tip



central open tip



dilatation tip



Couvelaire tip



Mercier tip



Dufour tip




Flute tip

Pictures may differ from the original



**2-way balloon catheter**

- ✓ In CH sizes 06 – 10 the catheter comes with a stabilising guiderail.
- ✓ Due to the high balloon capacity of up to 30 ml, the catheter is safely fixed in the bladder.
- ✓ Use our blocker syringe containing 10% glycerine solution here.
- ✓ Excellent biocompatibility and high resistance to encrustation.
- ✓ The X-ray line allows for radiological imaging of the catheter.
- ✓ MRT – compatibility.
- ✓ The catheter can be used for long-term drainage.

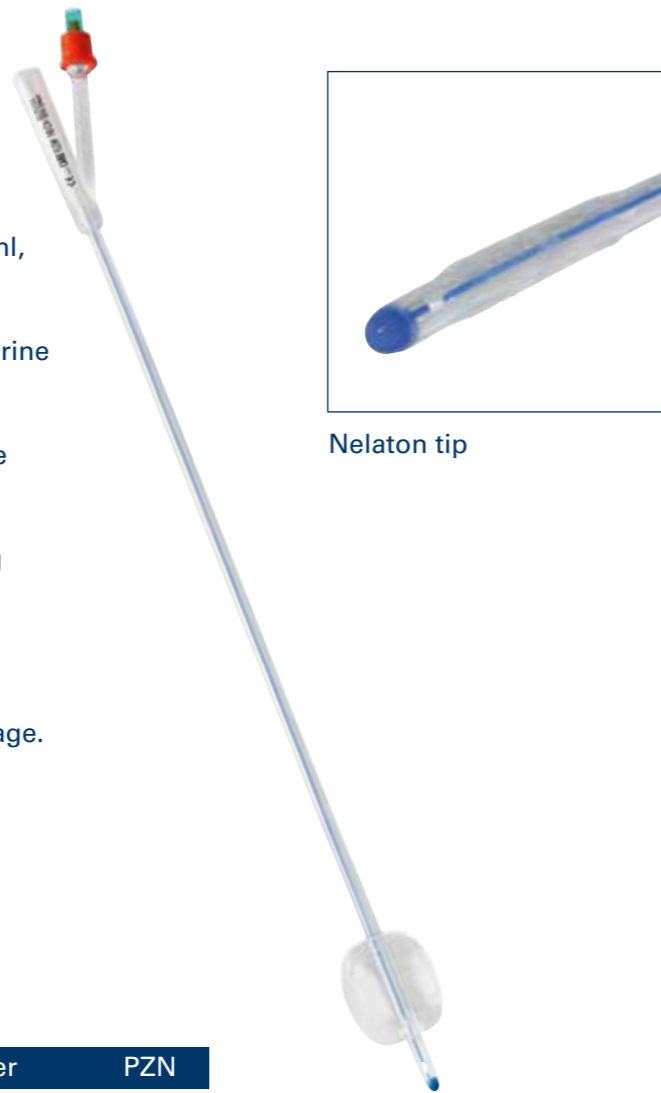

Nelaton tip



Balloon volume: 03-30 ml

Size	Balloon	Length	Article number	PZN
CH 06	03 ml	27 cm	01/01/11/06/03/GHC	01150196
CH 08	05 ml	27 cm	01/01/11/08/05/GHC	01150204
CH 10	05 ml	27 cm	01/01/11/10/05/GHC	01150210
CH 12	10 ml	40 cm	01/01/11/12/10/GHC	01150368
CH 14	10 ml	40 cm	01/01/11/14/10/GHC	01150374
CH 16	30 ml	40 cm	01/01/11/16/30/GHC	01150380
CH 18	30 ml	40 cm	01/01/11/18/30/GHC	01150397
CH 20	30 ml	40 cm	01/01/11/20/30/GHC	01150428
CH 22	30 ml	40 cm	01/01/11/22/30/GHC	01150434
CH 24	30 ml	40 cm	01/01/11/24/30/GHC	01150440
CH 26	30 ml	40 cm	01/01/11/26/30/GHC	01150457

**Packaging Unit:** 10 pieces  
**Packaging:** sterile

Pictures may differ from the original



**2-way balloon catheter**

- ✓ Very high wearing comfort, due to the soft-silicone material of the catheter shaft.
- ✓ Secure fixation of the catheter in the bladder due to the high balloon capacity of up to 30 ml
- ✓ Use our blocking syringe containing 10% glycerine solution here.
- ✓ Excellent biocompatibility and high resistance to encrustation.
- ✓ The x-ray line allows for radiological imaging of the catheter.
- ✓ The catheter can be used for long-term drainage.



Nelaton tip



Balloon volume: 03-30 ml

Size	Balloon	Length	Article number	PZN
CH 06	03 ml	31 cm	01/01/01/06/03/GHC	17883185
CH 08	05 ml	31 cm	01/01/01/08/05/GHC	17883216
CH 10	05 ml	31 cm	01/01/01/10/05/GHC	17883222
CH 12	10 ml	40 cm	01/01/01/12/10/GHC	17883239
CH 14	10 ml	40 cm	01/01/01/14/10/GHC	17883251
CH 16	30 ml	40 cm	01/01/01/16/30/GHC	17883268
CH 18	30 ml	40 cm	01/01/01/18/30/GHC	17883274
CH 20	30 ml	40 cm	01/01/01/20/30/GHC	17883280
CH 22	30 ml	40 cm	01/01/01/22/30/GHC	17883297
CH 24	30 ml	40 cm	01/01/01/24/30/GHC	17883305
CH 26	30 ml	40 cm	01/01/01/26/30/GHC	17883311


**Packaging Unit:** 10 pieces  
**Packaging:** sterile

Pictures may differ from the original



**2-way *integral* balloon catheter**

- ✓ In CH sizes 08-10, the catheters are equipped with a stabilizing guiderail.
- ✓ Secure fixation of the catheter in the bladder thanks to the *integral* balloon incorporated into the shaft.
- ✓ Use our blocker syringe containing 10% glycerine solution here.
- ✓ Excellent biocompatibility and high resistance to encrustation.
- ✓ The x-ray line allows for radiological imaging of the catheter.
- ✓ MRT – compatibility.
- ✓ The catheter can be used for long-term drainage.

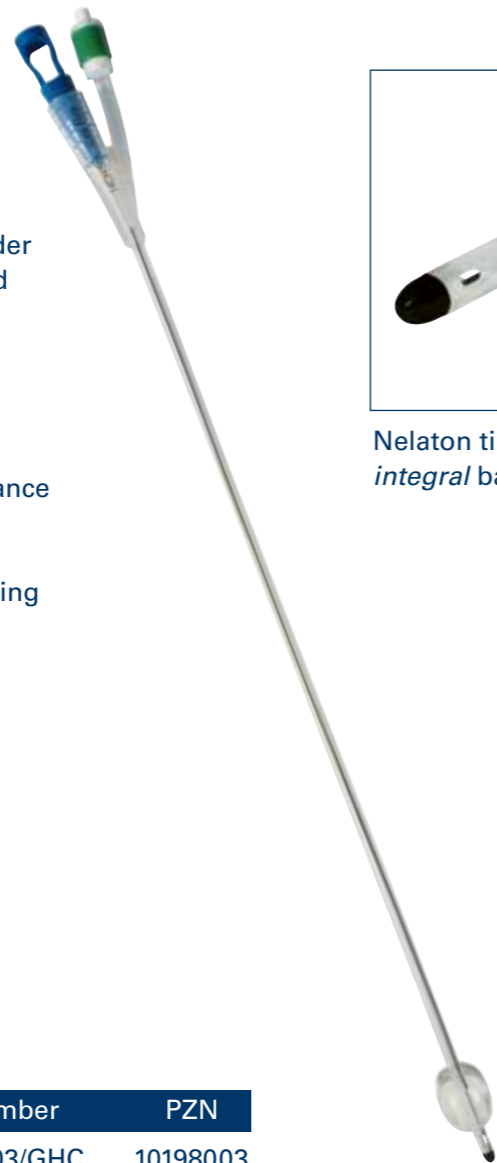

Nelaton tip, *integral* balloon



Balloon volume: 03-10 ml

Size	Balloon	Length	Article number	PZN
CH 08	03 ml	33 cm	01/01/23/08/03/GHC	10198003
CH 10	03 ml	33 cm	01/01/23/10/03/GHC	10198026
CH 12	05 ml	40 cm	01/01/23/12/05/GHC	10198032
CH 14	05 ml	40 cm	01/01/23/14/05/GHC	10198049
CH 16	10 ml	40 cm	01/01/23/16/10/GHC	10198055
CH 18	10 ml	40 cm	01/01/23/18/10/GHC	10198061
CH 20	10 ml	40 cm	01/01/23/20/10/GHC	10198078
CH 22	10 ml	40 cm	01/01/23/22/10/GHC	10198084
CH 24	10 ml	40 cm	01/01/23/24/10/GHC	10198090
CH 26	10 ml	40 cm	01/01/23/26/10/GHC	10198109


**Packaging Unit:** 10 pieces  
**Packaging:** sterile

Pictures may differ from the original



**2-way *integral* balloon catheter**

- ✓ Secure fixation of the catheter in the bladder thanks to the *integral* balloon incorporated into the shaft.
- ✓ Use our blocker syringe containing 10% glycerine solution here.
- ✓ Excellent biocompatibility and high resistance to encrustation.
- ✓ MRT – compatibility.
- ✓ The catheter can be used for long-term drainage.

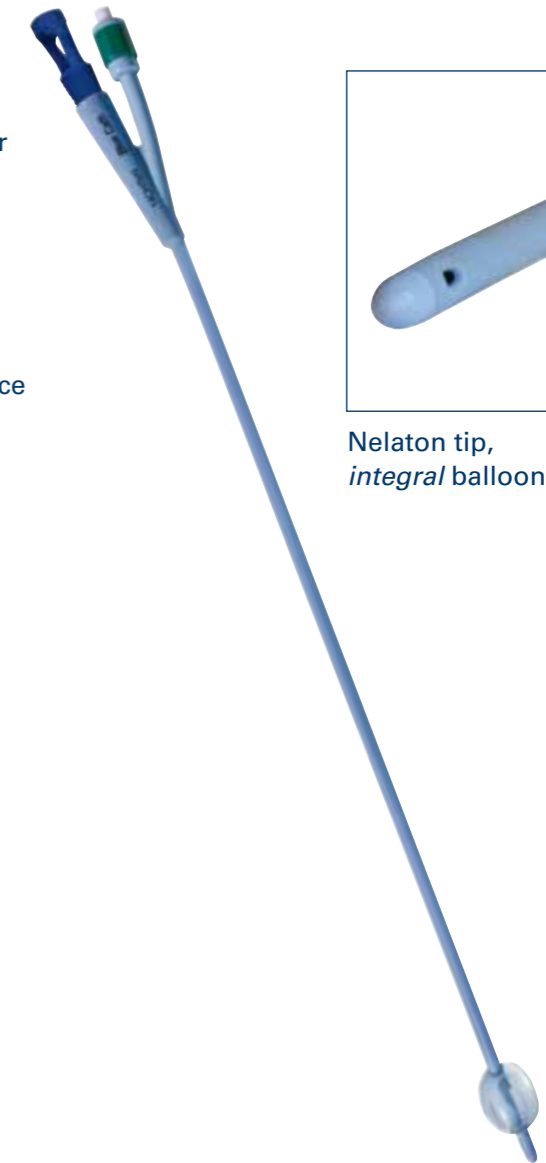

Nelaton tip, *integral* balloon



Balloon volume: 05-10 ml

Size	Balloon	Length	Article number	PZN
CH 12	05 ml	40 cm	01/01/24/12/05/GHC	11850871
CH 14	05 ml	40 cm	01/01/24/14/05/GHC	11850888
CH 16	10 ml	40 cm	01/01/24/16/10/GHC	11850894
CH 18	10 ml	40 cm	01/01/24/18/10/GHC	11850902
CH 20	10 ml	40 cm	01/01/24/20/10/GHC	11850919
CH 22	10 ml	40 cm	01/01/24/22/10/GHC	11850925
CH 24	10 ml	40 cm	01/01/24/24/10/GHC	11850931


**Packaging Unit:** 10 pieces  
**Packaging:** sterile

Pictures may differ from the original



2-way balloon catheter

- ✓ High elasticity and wearing comfort, due to the silicone coating on the inner and outer catheter wall.
- ✓ Due to the high balloon capacity, the catheter is safety fixed in the bladder.
- ✓ For blocking use our blocker syringe containing 10% glycerine solution.
- ✓ Excellent biocompatibility and high resistance to encrustation.
- ✓ The catheter can be used for long-term drainage.

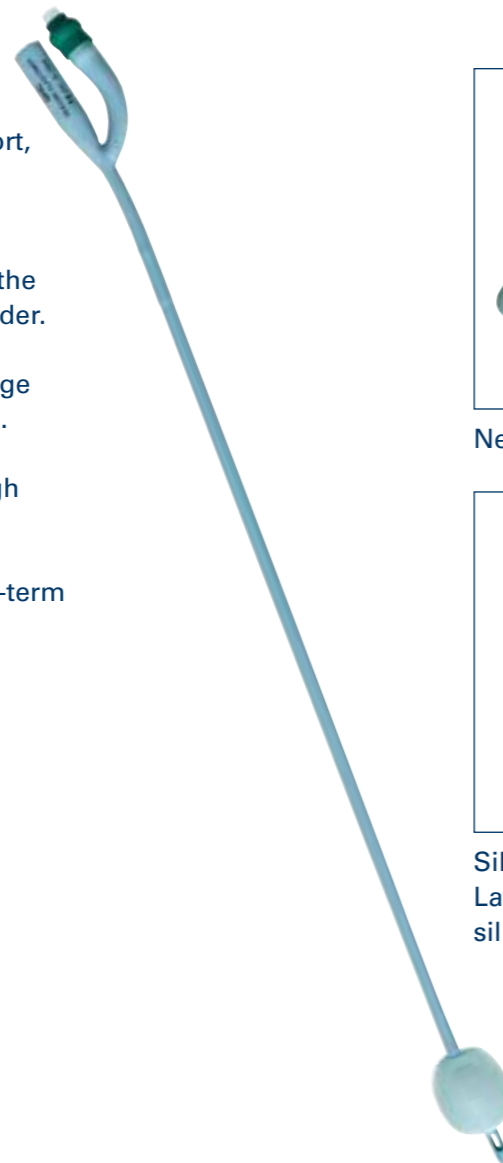

Nelaton tip



Silicone elastomer – Latex core with silicone coating



Balloon volume: 10-30 ml

Size	Balloon	Length	Article number	PZN
CH 12	10 ml	40 cm	01/01/10/12/10/GHC	01149371
CH 14	30 ml	40 cm	01/01/10/14/30/GHC	01149388
CH 16	30 ml	40 cm	01/01/10/16/30/GHC	01149394
CH 18	30 ml	40 cm	01/01/10/18/30/GHC	01149939
CH 20	30 ml	40 cm	01/01/10/20/30/GHC	01149945
CH 22	30 ml	40 cm	01/01/10/22/30/GHC	01149951
CH 24	30 ml	40 cm	01/01/10/24/30/GHC	01149968
CH 26	30 ml	40 cm	01/01/10/26/30/GHC	01150173


Packaging Unit: 10 pieces  
Packaging: sterile

Pictures may differ from the original



2-way balloon catheter

- ✓ An atraumatic dilating ball tip ensures a smooth passage of the prostatic fossa.
- ✓ Safe fixation of the catheter due to its high balloon capacity of up to 30 ml.
- ✓ Use our blocking syringe containing 10% glycerine solution here.
- ✓ Excellent biocompatibility and high resistance to encrustation.
- ✓ The x-ray line allows for radiological imaging of the catheter.
- ✓ MRT – compatibility.
- ✓ The catheter can be used for long term drainage.

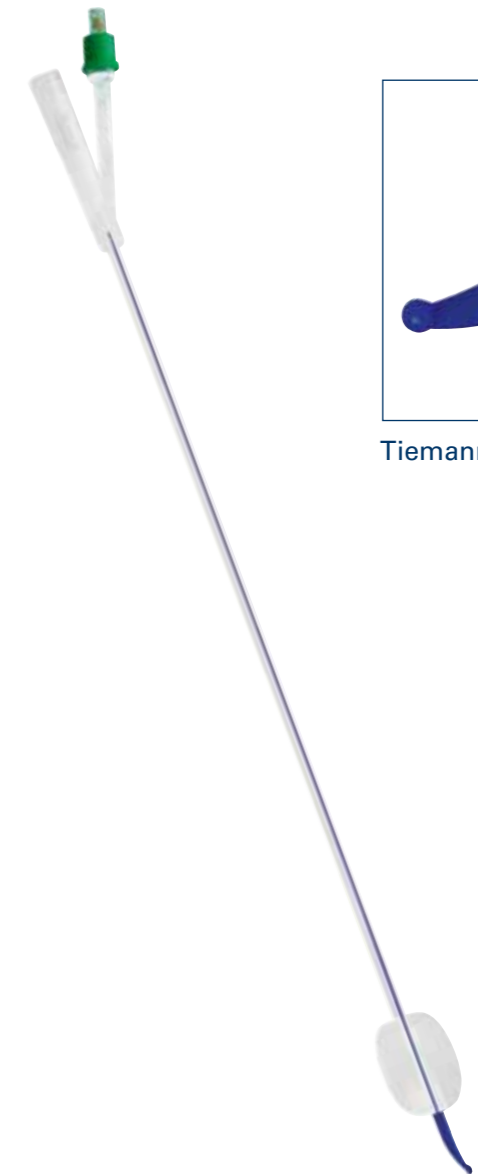

Tiemann tip



Balloon volume: 5-30 ml

Size	Balloon	Length	Article number	PZN
CH 12	05 ml	40 cm	01/01/08/12/05/GHC	17933670
CH 14	20 ml	40 cm	01/01/08/14/20/GHC	17933687
CH 16	30 ml	40 cm	01/01/08/16/30/GHC	17933693
CH 18	30 ml	40 cm	01/01/08/18/30/GHC	17933701
CH 20	30 ml	40 cm	01/01/08/20/30/GHC	17933718
CH 22	30 ml	40 cm	01/01/08/22/30/GHC	17933724


Packaging Unit: 10 pieces  
Packaging: sterile

Pictures may differ from the original



**2-way balloon catheter**

- ✓ Excellent drainage performance, due to the central open tip, drainage eye, and high lumen, even for problematic urine findings.
- ✓ An atraumatic flute tip ensures an easy placement of the catheter.
- ✓ Safe fixation of the catheter due to its high balloon capacity of up to 30 ml.
- ✓ Use our blocking syringe containing 10% glycerine solution here.
- ✓ Excellent biocompatibility and high resistance to encrustation.
- ✓ The x-ray line allows for radiological imaging of the catheter.
- ✓ MRT – compatibility.
- ✓ The catheter can be used for long term drainage.

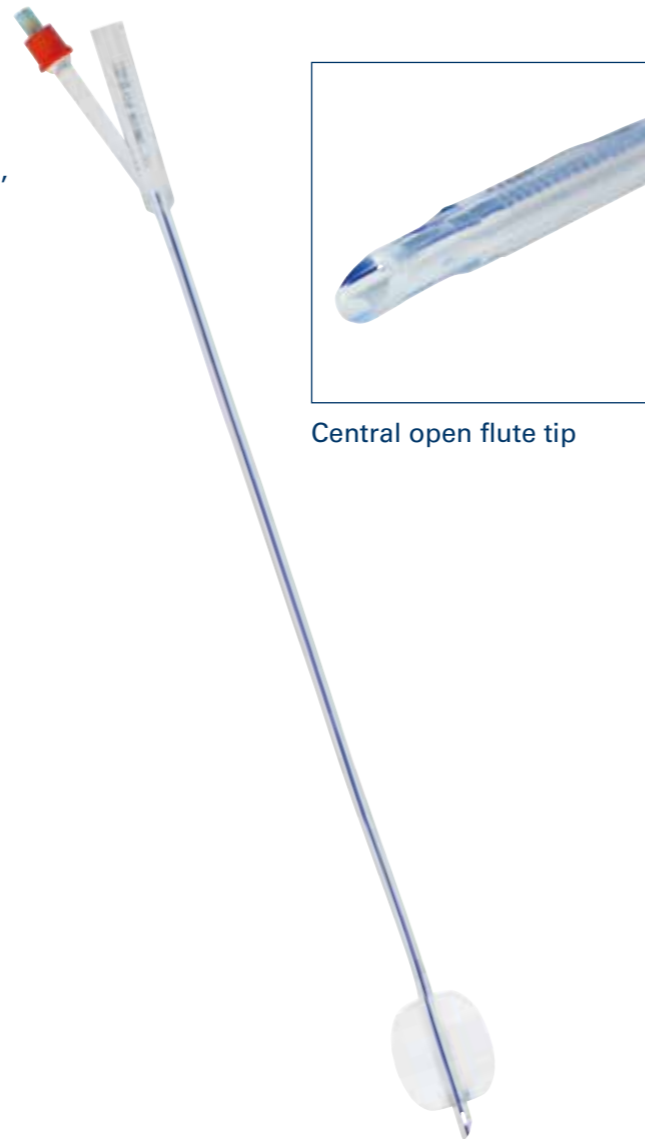

Central open flute tip



Balloon volume: 10-30 ml

Size	Balloon	Length	Article number	PZN
CH 12	10 ml	40 cm	01/01/07/12/10/GHC	17933747
CH 14	10 ml	40 cm	01/01/07/14/10/GHC	17933753
CH 16	30 ml	40 cm	01/01/07/16/30/GHC	17933776
CH 18	30 ml	40 cm	01/01/07/18/30/GHC	17933782
CH 20	30 ml	40 cm	01/01/07/20/30/GHC	17933799
CH 22	30 ml	40 cm	01/01/07/22/30/GHC	17933807



**Packaging Unit:** 10 pieces  
**Packaging:** sterile

Pictures may differ from the original



**2-way balloon catheter**

- ✓ High elasticity and wearing comfort, due to the silicone coating on the inner and outer catheter wall.
- ✓ Optimal drainage and irrigation due to two staggered drainage eyes and open couvelaire tip.
- ✓ Excellent biocompatibility and high resistance to encrustation.
- ✓ Safe fixation of the catheter, due to the high balloon capacity of up to 60 ml.
- ✓ For blocking use our blocker syringe containing 10% glycerine solution.
- ✓ The catheter can be used for long-term drainage.

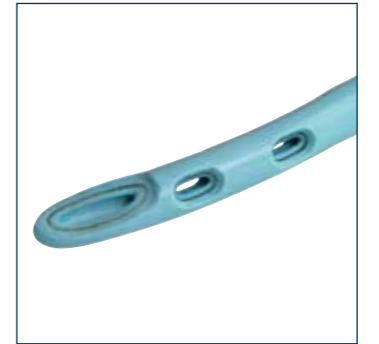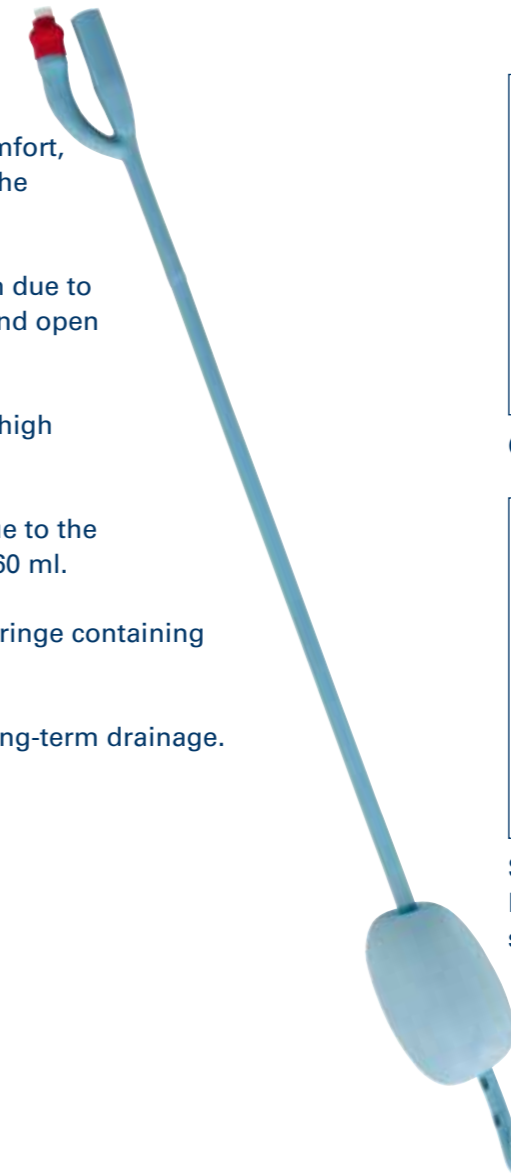

Couvelaire tip



Silicone elastomer – Latex core with silicone coating



Balloon volume: up to 60 ml

Size	Balloon	Length	Article number	PZN
CH 18	60 ml	40 cm	01/01/15/18/60/GHC	06434975
CH 20	60 ml	40 cm	01/01/15/20/60/GHC	06434952
CH 22	60 ml	40 cm	01/01/15/22/60/GHC	06434969
CH 24	60 ml	40 cm	01/01/15/24/60/GHC	09710080

**Packaging Unit:** 10 pieces  
**Packaging:** sterile

Pictures may differ from the original



**2-way balloon catheter**

- ✓ Made of high-quality, siliconized medical latex.
- ✓ The catheter offers high patient comfort.
- ✓ Due to the high balloon capacity, the catheter is safety fixed in the bladder.
- ✓ Use our blocker syringe containing 10% glycerine solution here.
- ✓ The catheter is suited for short-term drainage of up to 5 days only.
- ✓ This product contains natural latex which may cause allergic reactions.



Nelaton tip



Balloon volume: 10-30 ml

Size	Balloon	Length	Article number	PZN
CH 12	10 ml	40 cm	01/01/13/12/10/GHC	01150701
CH 14	10 ml	40 cm	01/01/13/14/10/GHC	01150859
CH 16	30 ml	40 cm	01/01/13/16/30/GHC	01150865
CH 18	30 ml	40 cm	01/01/13/18/30/GHC	01150871
CH 20	30 ml	40 cm	01/01/13/20/30/GHC	01150888
CH 22	30 ml	40 cm	01/01/13/22/30/GHC	01150919

**Packaging Unit:** 10 pieces  
**Packaging:** sterile

Pictures may differ from the original



**PVC single-use male catheter, transparent**

- ✓ Disposable male catheter with two staggered atraumatic drainage eyes.
- ✓ Rounded and closed Nelaton or Tiemann tip, for a smooth passage of the prostatic fossa.
- ✓ Soft PVC material.
- ✓ The catheter is for single use only.



Nelaton tip



Tiemann tip



Laterally offset eyes

**Nelaton tip**

Size	Length	Article number	PZN
CH 06	24 cm	01/01/35/06/00/GHC	18320969
CH 08	24 cm	01/01/35/08/00/GHC	18320975
CH 10	24 cm	01/01/35/10/00/GHC	18320981
CH 12	24 cm	01/01/35/12/00/GHC	18320998
CH 14	24 cm	01/01/35/14/00/GHC	18321006
CH 16	24 cm	01/01/35/16/00/GHC	18321035

Size	Length	Article number	PZN
CH 10	40 cm	01/01/36/10/00/GHC	18321058
CH 12	40 cm	01/01/36/12/00/GHC	18321064
CH 14	40 cm	01/01/36/14/00/GHC	18321070
CH 16	40 cm	01/01/36/16/00/GHC	18321087
CH 18	40 cm	01/01/36/18/00/GHC	18321093

**Tiemann tip**

Size	Length	Article number	PZN
CH 12	40 cm	01/01/37/12/00/GHC	18321124
CH 14	40 cm	01/01/37/14/00/GHC	18321130
CH 16	40 cm	01/01/37/16/00/GHC	18321147
CH 18	40 cm	01/01/37/18/00/GHC	18321153
CH 20	40 cm	01/01/37/20/00/GHC	18321176

**Packaging Unit:** 100 pieces  
**Packaging:** sterile

Pictures may differ from the original



**3-way balloon catheter**

- ✓ An atraumatic Nelaton tip ensures a smooth passage of the prostatic fossa.
- ✓ Optimal drainage and irrigation performance, due to two staggered drainage eyes.
- ✓ Excellent biocompatibility and high resistance to encrustation.
- ✓ Due to the high balloon capacity, the catheter is safely fixed in the bladder.
- ✓ Use our blocker syringe containing 10% glycerine solution here.
- ✓ The x-ray line allows for radiological imaging of the catheter.
- ✓ MRT – compatibility
- ✓ The catheter can be used for post-operative long-term drainage.



Nelaton tip



Balloon volume: 30 ml

Size	Balloon	Length	Article number	PZN
CH 18	30 ml	40 cm	01/01/22/18/30/GHC	15618237
CH 20	30 ml	40 cm	01/01/22/20/30/GHC	15618243
CH 22	30 ml	40 cm	01/01/22/22/30/GHC	15618266
CH 24	30 ml	40 cm	01/01/22/24/30/GHC	15618295

**Packaging Unit:** 10 pieces  
**Packaging:** sterile

Pictures may differ from the original



**3-way balloon catheter**

- ✓ An atraumatic and kink resistant Dufour tip ensures a smooth passage of the prostatic fossa.
- ✓ Optimal drainage and irrigation performance, due to two staggered drainage eyes and an open tip.
- ✓ Excellent biocompatibility and high resistance to encrustation.
- ✓ Due to the high balloon capacity, the catheter is safely fixed in the bladder.
- ✓ Use our blocker syringe containing 10% glycerine solution here.
- ✓ The x-ray tip allows for radiological imaging of the catheter.
- ✓ MRT – compatibility.
- ✓ The catheter can be used for post-operative long-term drainage.



Dufour tip



Balloon volume: 40-70 ml

Size	Balloon	Length	Article number	PZN
CH 18	40 ml	40 cm	01/01/16/18/40/GHC	01151178
CH 20	50 ml	40 cm	01/01/16/20/50/GHC	01151184
CH 22	60 ml	40 cm	01/01/16/22/60/GHC	01151190
CH 24	70 ml	40 cm	01/01/16/24/70/GHC	01151379

**Packaging Unit:** 10 pieces  
**Packaging:** sterile

Pictures may differ from the original



**3-way balloon catheter**

- ✓ Optimal drainage and irrigation due to two staggered drainage eyes and an open tip.
- ✓ An atraumatic and kink resistant Couvelaire tip ensures a smooth passage of the prostatic fossa.
- ✓ Excellent biocompatibility and high resistance to encrustation.
- ✓ Due to the high balloon capacity, the catheter is safely fixed in the bladder.
- ✓ Use our blocker syringe containing 10% glycerine solution here.
- ✓ The x-ray tip allows for radiological imaging of the catheter.
- ✓ MRT – compatibility.
- ✓ The catheter can be used for post-operative long-term drainage.



Couvelaire tip



Balloon volume: 40-70 ml

Size	Balloon	Length	Article number	PZN
CH 18	40 ml	40 cm	01/01/17/18/40/GHC	01151385
CH 20	50 ml	40 cm	01/01/17/20/50/GHC	01151391
CH 22	60 ml	40 cm	01/01/17/22/60/GHC	01151416
CH 24	70 ml	40 cm	01/01/17/24/70/GHC	01151439


**Packaging Unit:** 10 pieces  
**Packaging:** sterile

Pictures may differ from the original



**3-way balloon catheter**

- ✓ An atraumatic and kink resistant Mercier tip ensures a smooth passage of the prostatic fossa.
- ✓ Optimal drainage and irrigation due to two staggered drainage eyes and an open tip.
- ✓ Excellent biocompatibility and high resistance to encrustation.
- ✓ Due to the high balloon capacity, the catheter is safely fixed in the bladder.
- ✓ Use our blocker syringe containing 10% glycerine solution here.
- ✓ The x-ray tip allows for radiological imaging of the catheter.
- ✓ MRT – compatibility.
- ✓ The catheter can be used for post-operative long-term drainage.

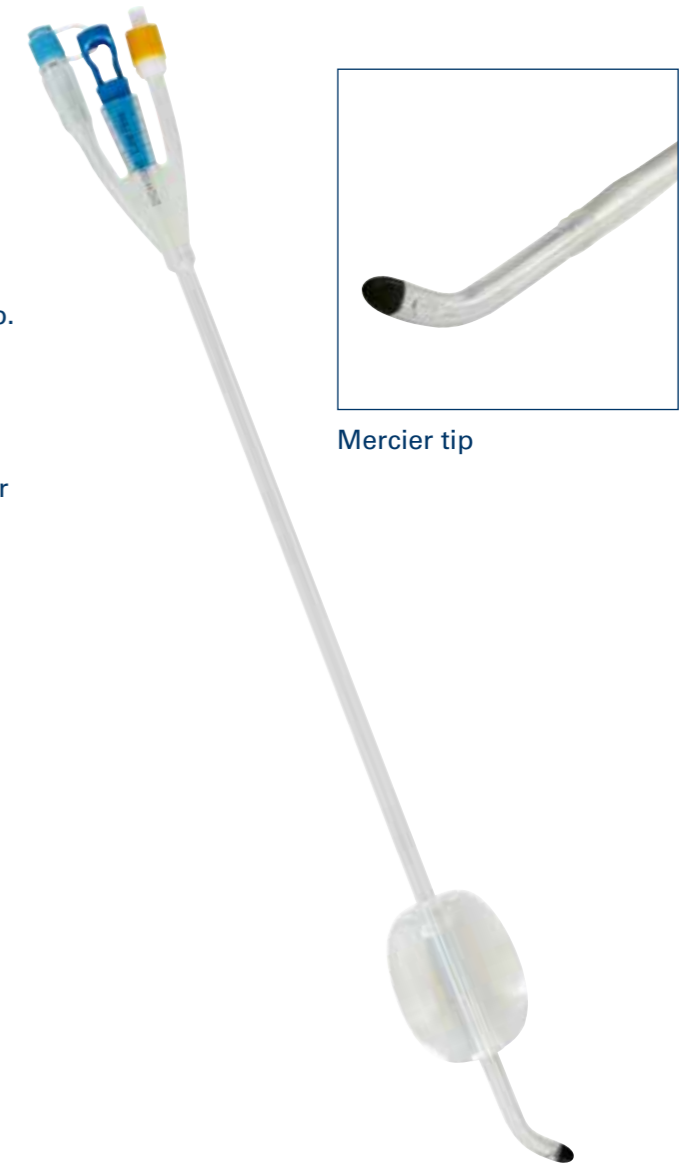

Mercier tip



Balloon volume: 40-70 ml

Size	Balloon	Length	Article number	PZN
CH 18	40 ml	40 cm	01/01/18/18/40/GHC	01151511
CH 20	50 ml	40 cm	01/01/18/20/50/GHC	01151528
CH 22	60 ml	40 cm	01/01/18/22/60/GHC	01151534
CH 24	70 ml	40 cm	01/01/18/24/70/GHC	01151540


**Packaging Unit:** 10 pieces  
**Packaging:** sterile

Pictures may differ from the original



**2-way balloon catheter**

- ✓ The catheter is equipped with an *integral* balloon, resulting in an even thickness throughout the catheter. This ensures that there is no painful compression of the fascia during catheter placement.
- ✓ The catheter offers high patient comfort.
- ✓ High drainage performance, due to the high lumen and the central open tip.
- ✓ We recommend using or blocking syringe containing a 10% glycerine solution, for a safe fixation of the catheter.
- ✓ Excellent biocompatibility and high resistance to encrustation.
- ✓ Short and central open tip.
- ✓ The x-ray line allows for radiological imaging of the catheter.
- ✓ MRT – compatibility.

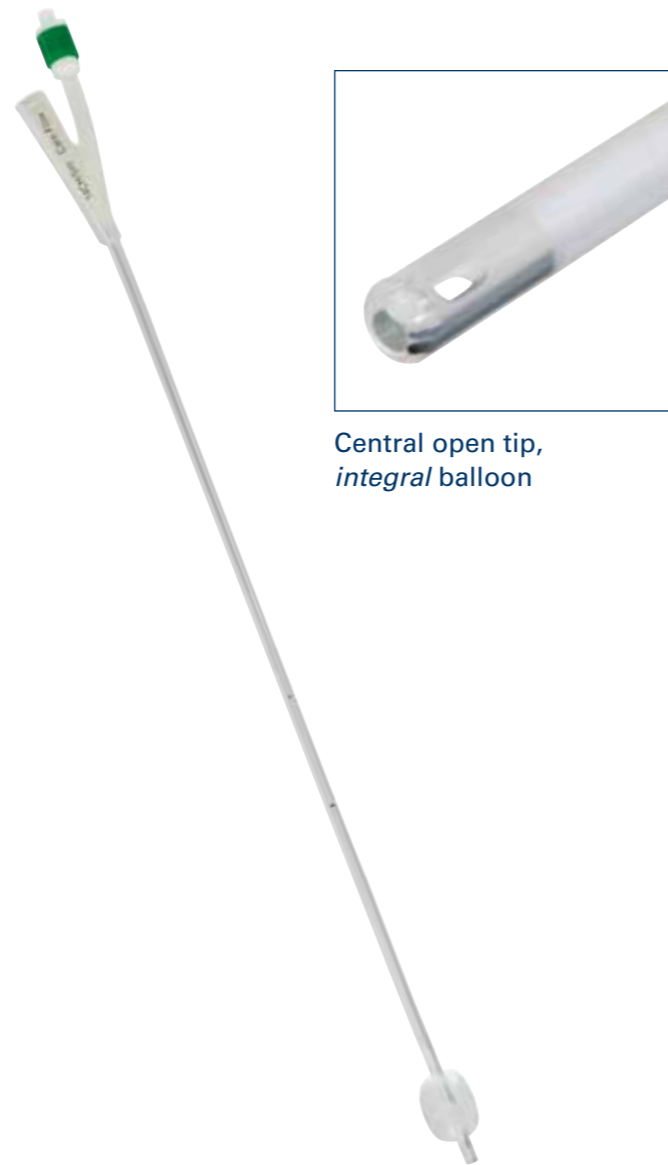

Central open tip,  
*integral* balloon



Balloonvolumen: 03-10 ml

Size	Balloon	Length	Article number	PZN
CH 10	03 ml	40 cm	01/02/12/10/03/GHC	01152568
CH 12	05 ml	40 cm	01/02/12/12/05/GHC	01152574
CH 14	05 ml	40 cm	01/02/12/14/05/GHC	01152605
CH 16	10 ml	40 cm	01/02/12/16/10/GHC	01152611
CH 18	10 ml	40 cm	01/02/12/18/10/GHC	01152628
CH 20	10 ml	40 cm	01/02/12/20/10/GHC	01152700
CH 22	10 ml	40 cm	01/02/12/22/10/GHC	01152717

**Packaging Unit:** 10 pieces  
**Packaging:** sterile

Pictures may differ from the original



**2-way *integral* balloon catheter**

- ✓ The catheter is equipped with an *integral* balloon, resulting in an even thickness throughout the catheter. This ensures that there is no painful compression of the fascia during catheter placement.
- ✓ The catheter offers high patient comfort.
- ✓ High drainage performance, due to the high lumen and the central open tip.
- ✓ We recommend using or blocking syringe containing a 10% glycerine solution, for a safe fixation of the catheter.
- ✓ Excellent biocompatibility and high resistance to encrustation.
- ✓ Short and central open tip.
- ✓ The x-ray line allows for radiological imaging of the catheter.
- ✓ MRT – compatibility.



Central open tip,  
*integral* balloon



Balloon volume: 03-10 ml

Size	Balloon	Length	Article number	PZN
CH 10	03 ml	40 cm	01/02/14/10/03/GHC	10198115
CH 12	05 ml	40 cm	01/02/14/12/05/GHC	10198121
CH 14	05 ml	40 cm	01/02/14/14/05/GHC	10198138
CH 16	10 ml	40 cm	01/02/14/16/10/GHC	10198144
CH 18	10 ml	40 cm	01/02/14/18/10/GHC	10198150
CH 20	10 ml	40 cm	01/02/14/20/10/GHC	10198167
CH 22	10 ml	40 cm	01/02/14/22/10/GHC	10198173


**Packaging Unit:** 10 pieces  
**Packaging:** sterile

Pictures may differ from the original



**2-way *integral* balloon catheter**

- ✓ The catheter is equipped with an *integral* balloon, resulting in an even thickness throughout the catheter. This ensures that there is no painful compression of the fascia during catheter placement.
- ✓ The catheter offers an excellent wearing comfort.
- ✓ We recommend using our blocking syringe containing a 10% glycerine solution, for a safe fixation of the catheter.
- ✓ 100% silicone for excellent biocompatibility and high resistance to encrustation.
- ✓ Short and central open tip.
- ✓ The x-ray line allows for radiological imaging of the catheter.
- ✓ MRT – compatibility.

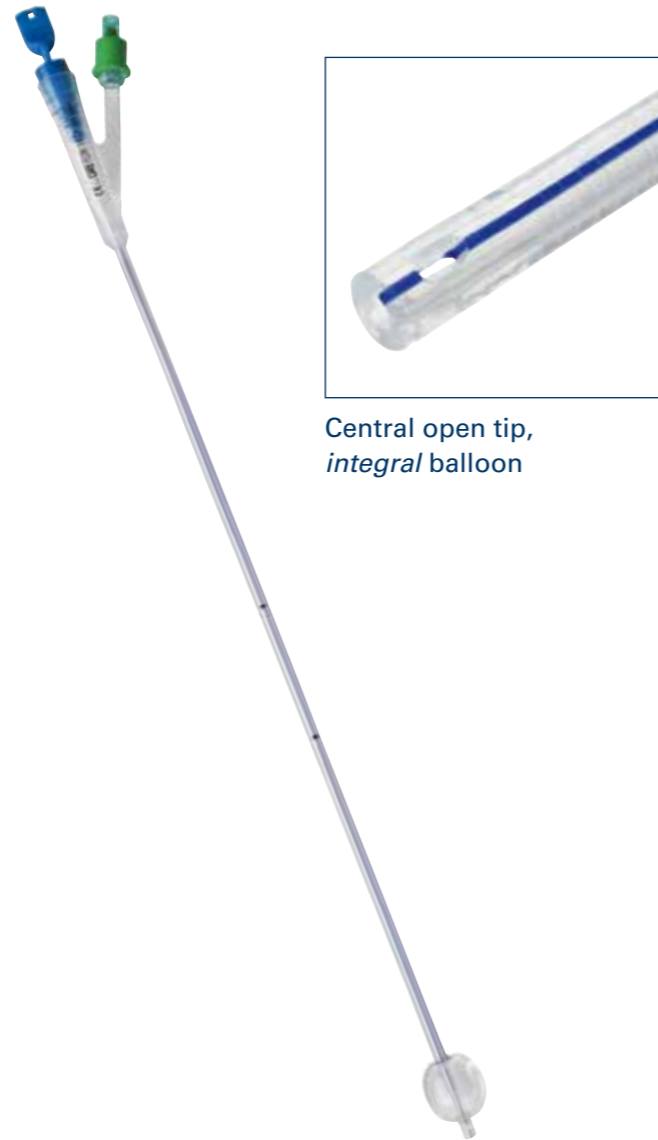

Central open tip, *integral* balloon



Balloon volume: 05-10 ml

Size	Balloon	Length	Article number	PZN
CH 12	05 ml	41 cm	01/02/09/12/05/GHC	16840937
CH 14	05 ml	41 cm	01/02/09/14/05/GHC	16840943
CH 16	10 ml	41 cm	01/02/09/16/10/GHC	16840972
CH 18	10 ml	41 cm	01/02/09/18/10/GHC	16840989
CH 20	10 ml	41 cm	01/02/09/20/10/GHC	16840995
CH 22	10 ml	41 cm	01/02/09/22/10/GHC	16841003

**Packaging Unit:** 10 pieces  
**Packaging:** sterile

Pictures may differ from the original



**2-way *integral* balloon catheter**

- ✓ The catheter is equipped with an *integral* balloon, resulting in an even thickness throughout the catheter. This ensures that there is no painful compression of the fascia during catheter placement.
- ✓ High drainage performance, due to the high lumen and the central open tip.
- ✓ Safe fixation of the catheter, due to the *integral* balloon.
- ✓ Use our blocker syringe containing 10% glycerine solution here.
- ✓ Excellent biocompatibility and high resistance to encrustation.
- ✓ Short and central open tip.
- ✓ The x-ray line allows for radiological imaging of the catheter.
- ✓ MRT – compatibility.



Central open tip, *integral* balloon



Balloon volume: 03-10 ml

Size	Balloon	Length	Article number	PZN
CH 10	03 ml	33 cm	01/02/13/10/03/GHC	01152746
CH 12	05 ml	33 cm	01/02/13/12/05/GHC	01152812
CH 14	05 ml	33 cm	01/02/13/14/05/GHC	01152829
CH 16	10 ml	33 cm	01/02/13/16/10/GHC	01152841
CH 18	10 ml	33 cm	01/02/13/18/10/GHC	01152918
CH 20	10 ml	33 cm	01/02/13/20/10/GHC	01152924
CH 22	10 ml	33 cm	01/02/13/22/10/GHC	01152947


**Packaging Unit:** 10 pieces  
**Packaging:** sterile

Pictures may differ from the original



**2-way *integral* balloon catheter**

- ✓ The catheter is equipped with an *integral* balloon, resulting in an even thickness throughout the catheter. This ensures that there is no painful compression of the fascia during catheter placement.
- ✓ Safe fixation of the catheter, due to the *integral* balloon.
- ✓ Use our blocker syringe containing 10% glycerine solution here.
- ✓ Excellent biocompatibility and high resistance to encrustation.
- ✓ Dilatation tip with central bore.
- ✓ MRT – compatibility.
- ✓ Graduation marks to ensure an optional placement.

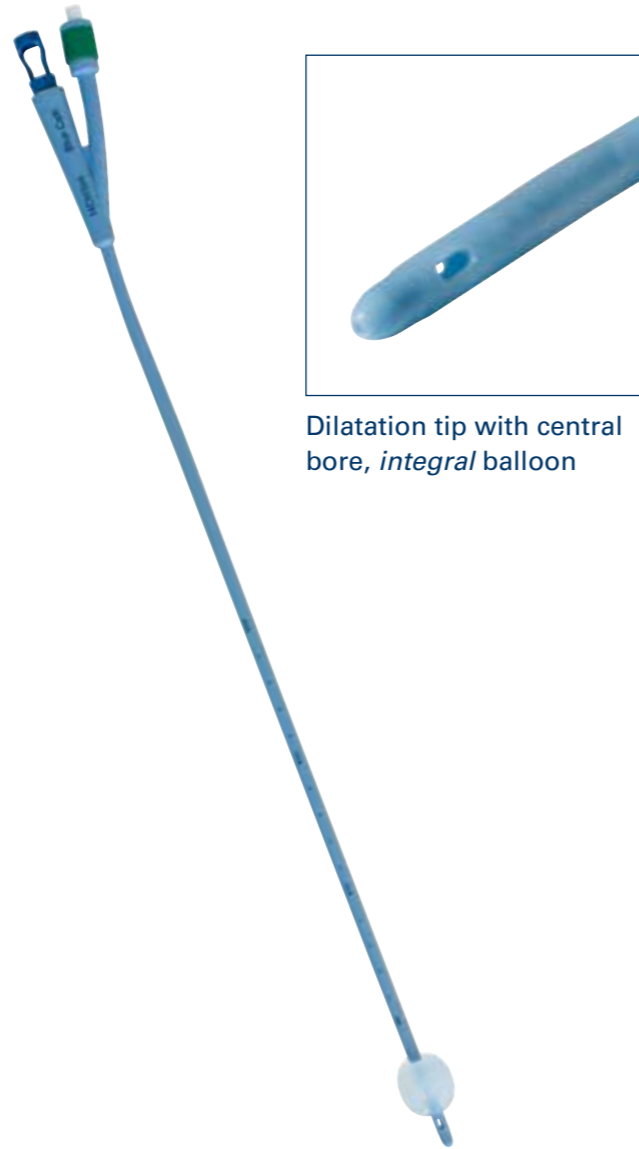

Dilatation tip with central bore, *integral* balloon



Balloon volume: 05-10 ml

Size	Balloon	Length	Article number	PZN
CH 12	05 ml	40 cm	01/02/21/12/05/GHC	11850948
CH 14	05 ml	40 cm	01/02/21/14/05/GHC	11850954
CH 16	10 ml	40 cm	01/02/21/16/10/GHC	11850983
CH 18	10 ml	40 cm	01/02/21/18/10/GHC	11851008
CH 20	10 ml	40 cm	01/02/21/20/10/GHC	11851014
CH 22	10 ml	40 cm	01/02/21/22/10/GHC	11851020


**Packaging Unit:** 10 pieces  
**Packaging:** sterile

Pictures may differ from the original



**2-way *integral* balloon catheter**

- ✓ The catheter is equipped with an *integral* balloon, resulting in an even thickness throughout the catheter. This ensures that there is no painful compression of the fascia during catheter placement.
- ✓ Safe fixation of the catheter, due to the *integral* balloon.
- ✓ High drainage performance, due to the high lumen and the central open tip.
- ✓ Use our blocker syringe containing 10% glycerine solution here.
- ✓ Excellent biocompatibility and high resistance to encrustation.
- ✓ Short, central open tip.
- ✓ MRT – compatibility.
- ✓ Graduation marks to ensure an optional placement.

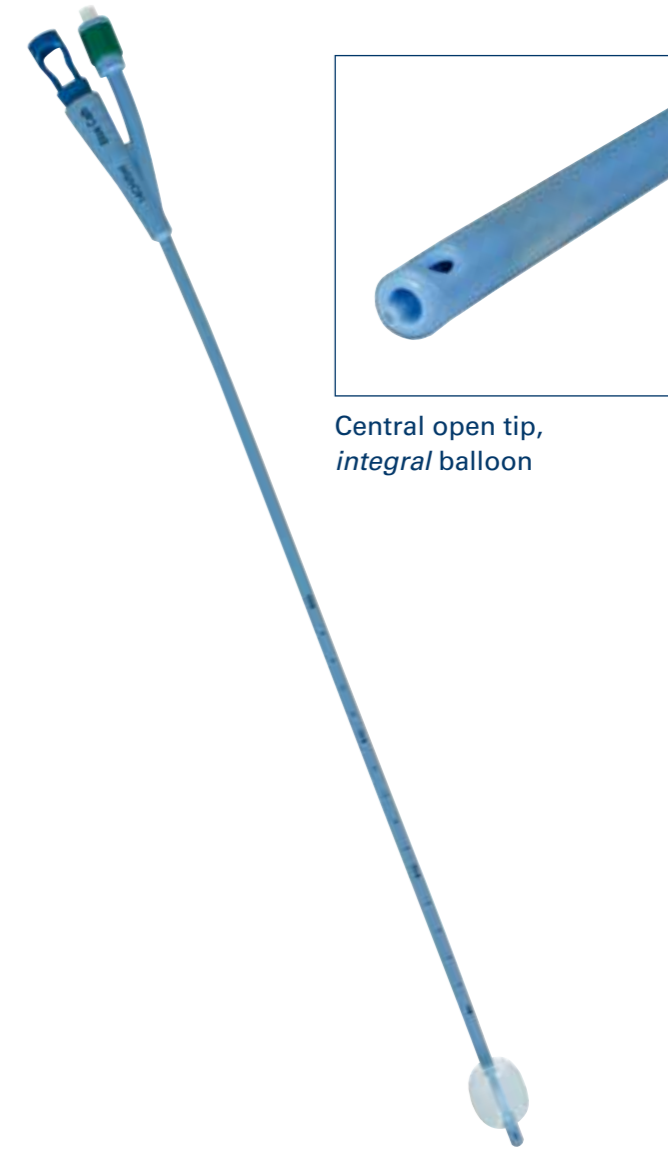

Central open tip, *integral* balloon



Balloon volume: 05-10 ml

Size	Balloon	Length	Article number	PZN
CH 12	05 ml	40 cm	01/02/22/12/05/GHC	11851037
CH 14	05 ml	40 cm	01/02/22/14/05/GHC	11851043
CH 16	10 ml	40 cm	01/02/22/16/10/GHC	11851072
CH 18	10 ml	40 cm	01/02/22/18/10/GHC	11851089
CH 20	10 ml	40 cm	01/02/22/20/10/GHC	11851095
CH 22	10 ml	40 cm	01/02/22/22/10/GHC	11851126
CH 24	10 ml	40 cm	01/02/22/24/10/GHC	11851132


**Packaging Unit:** 10 pieces  
**Packaging:** sterile

Pictures may differ from the original



**2-way *integral* blue balloon catheter**

- ✓ The catheter is equipped with an *integral* balloon, resulting in an even thickness throughout the catheter. This ensures that there is no painful compression of the fascia during catheter placement.
- ✓ The catheter offers an excellent wearing comfort.
- ✓ We recommend using our blocking syringe containing a 10% glycerine solution, for a safe fixation of the catheter.
- ✓ 100% silicone for excellent biocompatibility and high resistance to encrustation.
- ✓ Short and central open tip.
- ✓ The x-ray line allows for radiological imaging of the catheter.
- ✓ MRT – compatibility.
- ✓ Graduation marks to ensure an optimal placement.

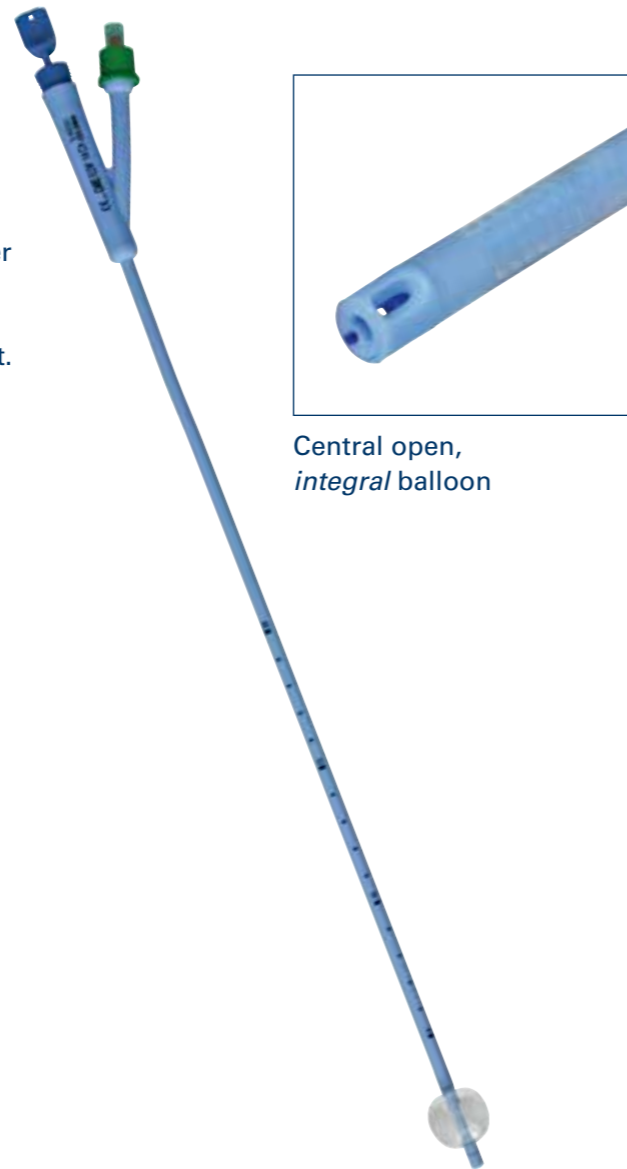

Central open, *integral* balloon



Balloon volume: 05-10 ml

Size	Balloon	Length	Article number	PZN
CH 12	05 ml	41 cm	01/02/11/12/05/GHC	16841026
CH 14	05 ml	41 cm	01/02/11/14/05/GHC	16841032
CH 16	10 ml	41 cm	01/02/11/16/10/GHC	16841049
CH 18	10 ml	41 cm	01/02/11/18/10/GHC	16841061
CH 20	10 ml	41 cm	01/02/11/20/10/GHC	16841078
CH 22	10 ml	41 cm	01/02/11/22/10/GHC	16841084


**Packaging Unit:** 10 pieces  
**Packaging:** sterile

Pictures may differ from the original



**2-way *integral* balloon catheter**

- ✓ The catheter is equipped with an *integral* balloon, resulting in an even thickness throughout the catheter. This ensures that there is no painful compression of the fascia during catheter placement.
- ✓ High drainage performance, due to the high lumen and the central open tip.
- ✓ Safe fixation of the catheter, due to the *integral* balloon.
- ✓ Use our blocker syringe containing 10% glycerine solution here.
- ✓ Excellent biocompatibility and high resistance to encrustation.
- ✓ Short, central open tip.
- ✓ MRT – compatibility.
- ✓ Graduation marks to ensure an optional placement.

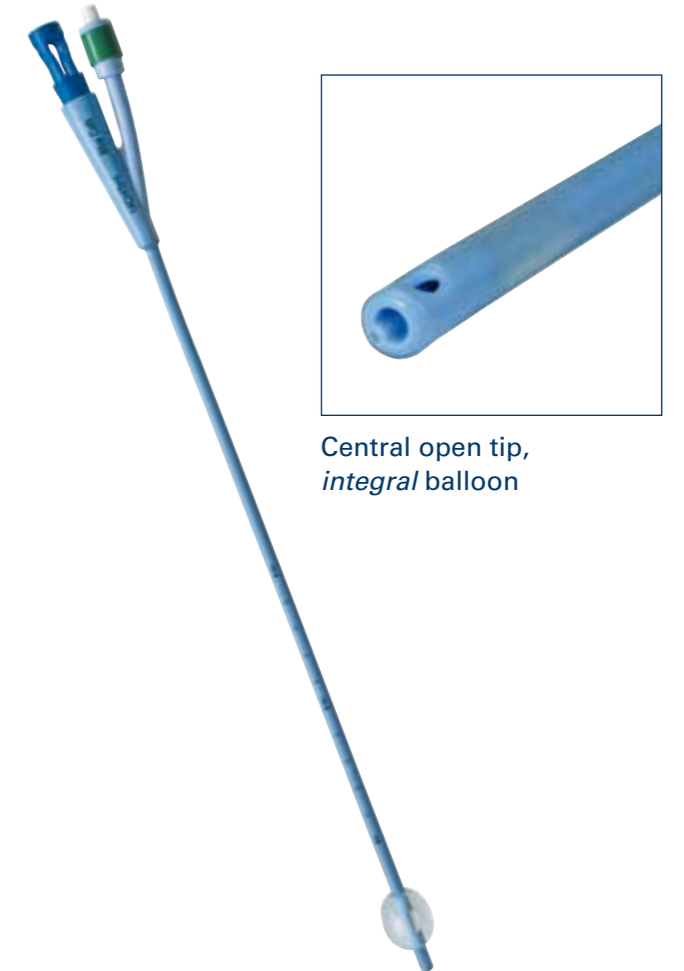

Central open tip, *integral* balloon



Balloon volume: 05-10 ml

Size	Balloon	Length	Article number	PZN
CH 12	05 ml	33 cm	01/02/23/12/05/GHC	11851149
CH 14	05 ml	33 cm	01/02/23/14/05/GHC	11851155
CH 16	10 ml	33 cm	01/02/23/16/10/GHC	11851161
CH 18	10 ml	33 cm	01/02/23/18/10/GHC	11851178
CH 20	10 ml	33 cm	01/02/23/20/10/GHC	11851184
CH 22	10 ml	33 cm	01/02/23/22/10/GHC	11851190

**Packaging Unit:** 10 pieces  
**Packaging:** sterile

Pictures may differ from the original




**2-way integral balloon catheter**

The puncture set contains:

- 100% silicone *integral* balloon catheter
- Stainless steel cannula with special cut
- Scalpel with small handle

- ✓ *Integral* balloon catheter with seamlessly incorporated balloon into the shaft.
- ✓ Optimal use of the inside diameter of the cannula due to *integral* balloon catheter: Charrière 10 and 12 with split cannula CH 13.8 (4.6 mm), puncture length 120 mm. Charrière 14 with split cannula CH 16.8 (5.6 mm), puncture length 120 mm.
- ✓ High drainage performance, due to the high lumen and the central open tip.
- ✓ Safe fixation of the catheter, due to the *integral* balloon.
- ✓ Use our blocker syringe containing 10% glycerine solution here.
- ✓ The x-ray line allows for radiological imaging of the catheter.
- ✓ MRT – compatibility.

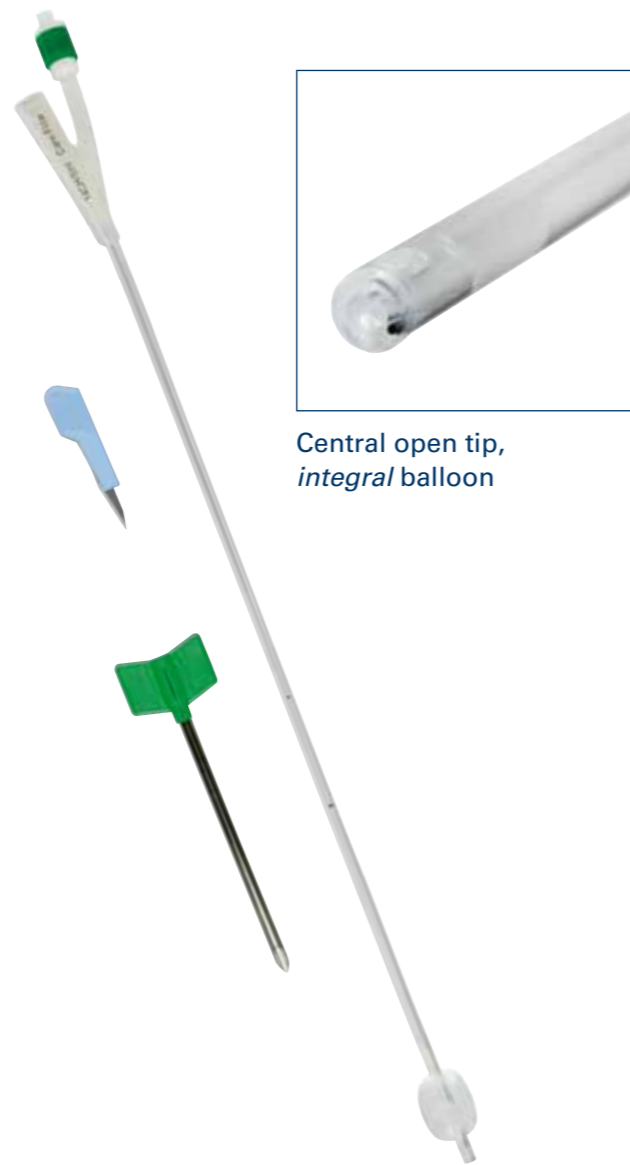

Central open tip, *integral* balloon



Balloon volume: 03-05 ml

Size	Balloon	Length	Article number	PZN
CH 10	03 ml	40 cm	01/02/15/10/03/GHC	01153059
CH 12	05 ml	40 cm	01/02/15/12/05/GHC	01153065
CH 14	05 ml	40 cm	01/02/15/14/05/GHC	01153088

**Packaging Unit:** 5 pieces  
**Packaging:** sterile




**2-way integral balloon catheter**

This charge set contains:

- 100% silicone *integral* balloon catheter
- Teflon-coated guidewire, 80 cm

- ✓ *Integral* balloon catheter with seamlessly incorporated balloon into the shaft.
- ✓ High drainage performance, due to the high lumen and the central open tip.
- ✓ Stiffened, Teflon-coated exchange guidewire, 80 cm in length, with a 3 cm straight, atraumatic flexible tip; wire diameter 0,035"
- ✓ Safe fixation of the catheter, due to the *integral* balloon.
- ✓ Use our blocker syringe containing 10% glycerine solution here.
- ✓ The x-ray line allows for radiological imaging of the catheter.
- ✓ MRT – compatibility.

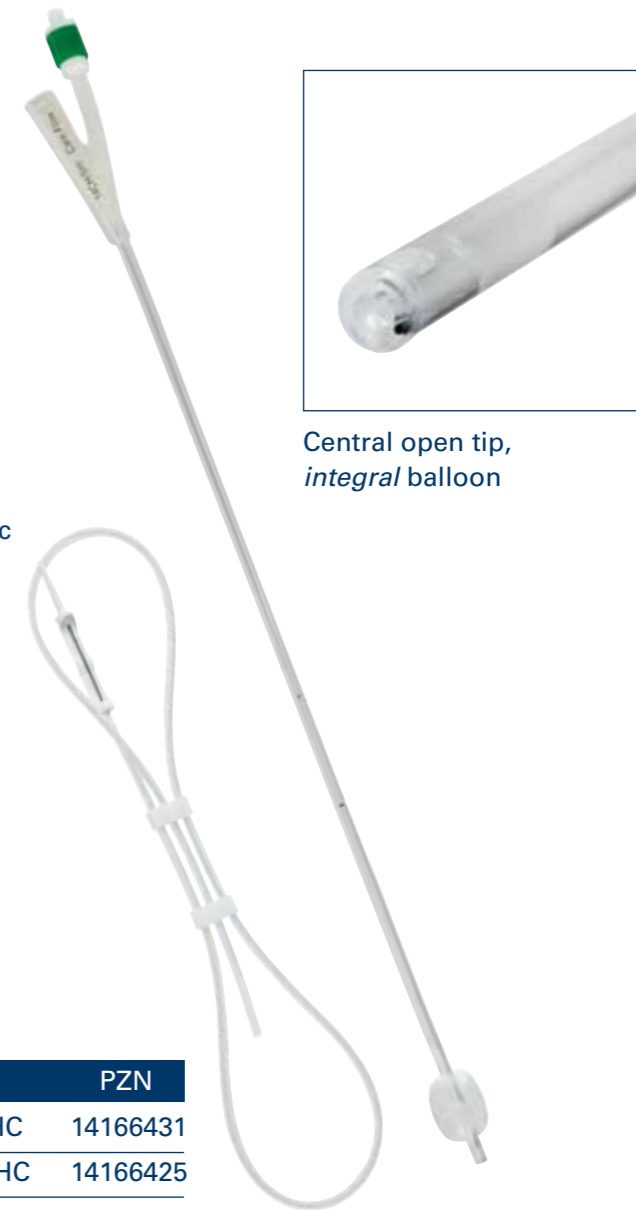

Central open tip, *integral* balloon



Balloon volume: 03-10 ml

**Guidewire**

Wire	Tip	Length	Article number	PZN
coated	3 cm flex.	80 cm	01/02/18/80/35/GHC	14166431
coated	3 cm flex.	100 cm	01/02/18/100/35/GHC	14166425

**Packaging Unit:** 1 pieces  
**Packaging:** sterile

**Change set, incl. Guidewire**

Size	Balloon	Length	Article number	PZN
CH 10	03 ml	40 cm	01/02/16/10/03/GHC	01153154
CH 12	05 ml	40 cm	01/02/16/12/05/GHC	01153160
CH 14	05 ml	40 cm	01/02/16/14/05/GHC	01153183
CH 16	10 ml	40 cm	01/02/16/16/10/GHC	16506380

**Packaging Unit:** 10 pieces  
**Packaging:** sterile



**2-way *integral* balloon catheter**

- ✓ Excellent biocompatibility and high resistance to encrustation.
- ✓ The catheters are suitable for long-term drainage.
- ✓ Dilating nelaton tip with central bore.
- ✓ High drainage performance, due to four drainage eyes located in front of and behind the balloon.
- ✓ The x-ray line allows for radiological imaging of the catheter.
- ✓ MRT – compatibility.
- ✓ Graduation marks to ensure an optimal placement.
- ✓ 100% silicone.



Dilating nelaton tip with central bore, *integral* balloon



Balloon volume: 1,5 - 5,0 ml

Size	Balloon	Length	Article number	PZN
CH 10	1,5 ml	40 cm	01/03/10/10/01/GHC	13851171
CH 12	2,0 ml	40 cm	01/03/10/12/02/GHC	13851188
CH 14	3,0 ml	40 cm	01/03/10/14/03/GHC	13851202
CH 16	5,0 ml	40 cm	01/03/10/16/05/GHC	13851219
CH 18	5,0 ml	40 cm	01/03/10/18/05/GHC	13851225

**Packaging Unit:** 10 pieces  
**Packaging:** sterile

Pictures may differ from the original



**Catheter insertion kit**

- ✓ The Basic Set comprises an inner and an outer package.
- ✓ The outer package of the basic set contains:  
one sterile 10 ml syringe for unblocking the placed catheter.
- ✓ The items in the inner, sterilised package are used to place a new transurethral catheter.



Article	Article number	PZN
Basic	01/01/99/BS/GHC	01972188
Clinic	01/01/99/BSK/GHC	08845872

**Components of the catheter insertion kit are:**

1 Tray 	6 Swabs 	1 Pair of tweezers 	2 Underpads 
1 Blocker syringe with 10% glycerine 	3 Nitrile gloves 	1 Fenestrated drape 	

**Packaging Unit:** 20 pieces  
**Packaging:** sterile

Additional included in Basic set:



10 ml syringe

Pictures may differ from the original



**Catheter insertion kit**

- ✓ The Catheter insertion Kits comprises an inner and an outer package.
- ✓ The outer package of the basic set contains: one sterile 10 ml syringe for unblocking the inserted catheter a mucous membrane disinfectant Octenisept 15 ml.
- ✓ Additional included in Comfort Plus: lubricant Endosgel 06 ml.
- ✓ The items in the inner, sterilised package are used to place a new transurethral catheter.



Article	Article number	PZN
Comfort	01/01/99/KL/CO/GHC	14209815
Comfort Plus	01/01/99/KL/CP/GHC	14209821

**Components of the catheter insertion kit are:**

1 Tray	6 Swabs	1 Pair of tweezers	10 ml Syringe	Mucous membrane disinfectant	Additional included in Comfort Plus:
					
1 Blocker syringe with 10% glycerine	3 Nitrile gloves	1 Fenestrated drape	2 Underpads	Octenisept 15 ml	Lubricant Endosgel 06 ml
					

**Packaging Unit:** 35 pieces  
**Packaging:** sterile

Pictures may differ from the original



**Catheter insertion kit**

- ✓ The Comfort Eco Set comprises an inner and an outer package.
- ✓ The outer package of the basic set contains: one sterile 10 ml syringe for unblocking the inserted catheter, mucous membrane disinfectant Octenisept 15 ml, lubricant LUBRICARE 6 ml.
- ✓ The items in the inner, sterilised package are used to place a new transurethral catheter.



Article	Article number	PZN
Comfort Eco	01/01/99/KL/CE/GHC	16134211

**Components of the catheter insertion kit are:**

1 Tray	6 Swabs	1 Pair of tweezers	10 ml Syringe	Mucous membrane disinfectant	Lubricant
					
1 Blocker syringe with 10% glycerine	3 NITRILE gloves	1 Fenestrated drape	2 Underpads	Octenisept, 15 ml	LUBRICARE, 06 ml
					

**Packaging Unit:** 35 pieces  
**Packaging:** sterile

Pictures may differ from the original




### Accessories for the catheter replacement

In addition to our catheters, drainage systems, and support devices, we supply all the necessary accessories to optimize continence care and catheter changes, ensuring exceptional service quality. Highlights include:

We want to highlight:

- **10 ml Block-Syringe with 10 % Glycerin Solution** (Art.-No. 01/01/99/KB/10/GHC): Maintains balloon volume throughout the indwelling period without diffusion of the fill medium. This keeps the balloon catheter securely in place, virtually eliminating the risk of dislocation.
- **LUBRICARE Lubricating Gel, 6 ml** (Art.-No. 01/01/99/GG/06/001): Adheres exceptionally well to mucosal surfaces and provides excellent lubrication, ensuring a smooth, atraumatic catheter exchange.

The products listed represent a selection from our range; additional items are available upon request.

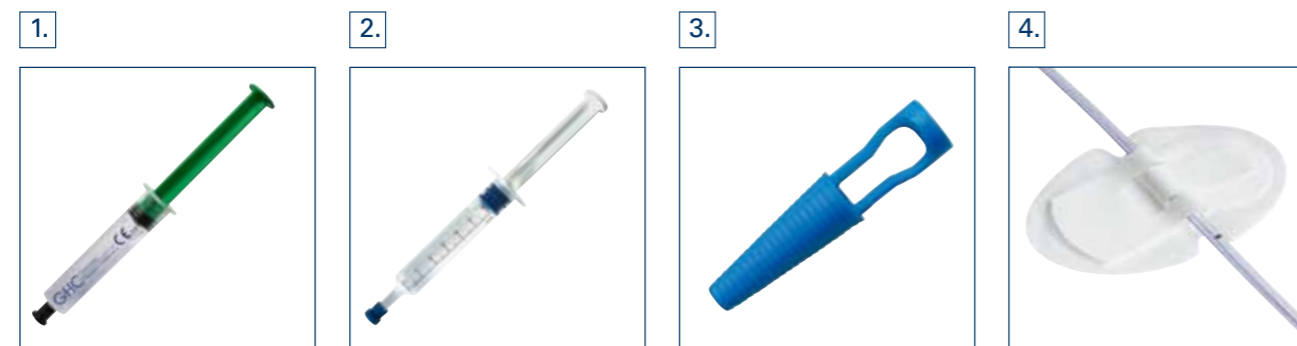

1. Blocker syringe with 10% glycerine, 10 ml

2. Lubricant LUBRICARE, 06 ml

3. Catheter plug

4. Velcro fixation plaster

Article	Article number	PZN	VE
1. Blocker syringe	01/01/99/KB/10/GHC	06718980	1VE = 200 pieces
2. LUBRICARE Lubricant	01/01/99/GG/06/001	13360220	1VE = 100 pieces
LUBRICARE Lubricant	01/01/99/GG/06/25	16938435	1VE = 4 x 25 pieces
3. Catheter plug	01/01/99/ZB/ST/GHC	14354919	1VE = 20 pieces
4. Velcro fixation plaster	08/02/FP/GHC	19511660	1VE = 10 pieces

Verpackung: steril

Pictures may differ from the original



### CARE FLOW catheter valve

- ✓ The optimal solution for mobile, independent users wearing a transurethral or suprapubic catheter.
- ✓ The catheter valve continuously trains the bladder muscles so that the bladder function remains completely intact. The emptying of the urinary bladder follows the user's personal rhythm. This sustainably enhances the quality of life.
- ✓ During the night the valve can be connected to a urine bag system by means of a firmly attached silicone tube.
- ✓ By flipping the lever of the flip valve by 180° the CARE FLOW catheter valve is completely opened as well as fixed in a safe position. This prevents an unintentional closing of the valve during the night.



Arrows indicate open / closed position

Closed catheter valve

Open catheter valve, flipped by 180°

Weight	Article number	PZN
3g	03/02/10/00/GHC	12378382

Packaging Unit: 5 pieces

Packaging: sterile

### CARE VENT catheter valve

- ✓ The special design of the CARE VENT prevents urine from dripping after the valve closes automatically. This ensures hygienic handling.
- ✓ The color coding of the rotary and slide controls makes one-handed movement extremely easy.



Weight	Article number	PZN
6g	03/02/10/01/GHC	01914676

Packaging Unit: 10 pieces

Packaging: sterile

Pictures may differ from the original

**CARE CON comfort | Typ I**  
**Male External Catheter (MEC) 100% silicone**



**Urinal condom, Self-adhesive**

- ✓ The CARE CON comfort captivates through its comfortable silicone material and good hold.
- ✓ The material thickness of the urinary sheath has been developed to cling softly and flexibly whilst ensuring a secure hold, which provides high safety.
- ✓ The stable connecting tube prevents the MEC from bending and ensures a secure connection to the urine bag.
- ✓ Small steps between the sizes ensures a perfect fit.
- ✓ The medical grade silicone maintains the breathing activity of the skin.
- ✓ Due to the use of transparent materials, skin alterations or irritations can be observed at any time.
- ✓ With a shaft length of 8,0 cm and an adhesive area of 4,5 cm, the CARE CON comfort meets the requirements of most MEC-users.



**CARE CON special | Typ II**  
**Male External Catheter (MEC) 100% silicone**



**Urinal condom, Self-adhesive**

- ✓ The CARE CON special captivates through its comfortable silicone material and good hold, despite the short shaft length.
- ✓ The material thickness of the urinary sheath has been developed to cling softly and flexibly whilst ensuring a secure hold, which provides high safety.
- ✓ The stable connecting tube prevents the MEC from bending and ensures a secure connection to the urine bag.
- ✓ Small steps between the sizes ensures a perfect fit.
- ✓ With a shaft length of 5,0 cm and an adhesive area of 4,0 cm, the CARE CON special meets the special requirements of users with a retracted penis.
- ✓ The medical grade silicone maintains the breathing activity of the skin.
- ✓ Due to the use of transparent materials, skin alterations or irritations can be observed at any time.



Diameter	Shaft / Adhesive surface	Article number	PZN
25 mm	8,0 cm / 4,5 cm	01/01/28/01/25/GHC	17520554
28 mm	8,0 cm / 4,5 cm	01/01/28/01/28/GHC	17520583
30 mm	8,0 cm / 4,5 cm	01/01/28/01/30/GHC	17520620
32 mm	8,0 cm / 4,5 cm	01/01/28/01/32/GHC	17520637
36 mm	8,0 cm / 4,5 cm	01/01/28/01/36/GHC	17520643
41 mm	8,0 cm / 4,5 cm	01/01/28/01/41/GHC	17520666

**Packaging Unit:** 30 pieces  
**Packaging:** non-sterile

Pictures may differ from the original

Diameter	Shaft / Adhesive surface	Article number	PZN
25 mm	5,0 cm / 4,0 cm	01/01/28/02/25/GHC	17520672
28 mm	5,0 cm / 4,0 cm	01/01/28/02/28/GHC	17520689
30 mm	5,0 cm / 4,0 cm	01/01/28/02/30/GHC	17520703
32 mm	5,0 cm / 4,0 cm	01/01/28/02/32/GHC	17520726
36 mm	5,0 cm / 4,0 cm	01/01/28/02/36/GHC	17520732
41 mm	5,0 cm / 4,0 cm	01/01/28/02/41/GHC	17520749

**Packaging Unit:** 30 pieces  
**Packaging:** non-sterile

Pictures may differ from the original

## CARE CON active | Typ III Male External Catheter (MEC) 100% silicone



### Urinal condom, Self-adhesive

- ✓ The CARE CON active captivates through its comfortable silicone material and good hold, with an extra-long adhesive area.
- ✓ The material thickness of the urinary sheath has been developed to cling softly and flexibly whilst ensuring a secure hold, which provides high safety.
- ✓ The stable connecting tube prevents the MEC from bending and ensures a secure connection to the urine bag.
- ✓ Small steps between the sizes ensures a perfect fit.
- ✓ With a shaft length of 8,0 cm and an adhesive area of 7,0 cm, the CARE CON active meets the requirements of sporty and active MEC-users.
- ✓ The medical grade silicone maintains the breathing activity of the skin.
- ✓ Due to the use of transparent materials, skin alterations or irritations can be observed at any time.



Diameter	Shaft / Adhesive surface	Article number	PZN
25 mm	8,0 cm / 7,0 cm	01/01/28/03/25/GHC	17520755
28 mm	8,0 cm / 7,0 cm	01/01/28/03/28/GHC	17520761
30 mm	8,0 cm / 7,0 cm	01/01/28/03/30/GHC	17520784
32 mm	8,0 cm / 7,0 cm	01/01/28/03/32/GHC	17520790
36 mm	8,0 cm / 7,0 cm	01/01/28/03/36/GHC	17520809
41 mm	8,0 cm / 7,0 cm	01/01/28/03/41/GHC	17520815

**Packaging Unit:** 30 pieces  
**Packaging:** non-sterile

Pictures may differ from the original

## CARE FLOW UB 4000TUR Closed irrigation bag system with pump ball for controlled post-operative bladder irrigation



### Irrigation bag system, 4.000 ml

- ✓ Can be fixed safely and correctly to many different types of bed frames by means of the twin holder with fixing cord.
- ✓ Bacteria- and fluid-tight venting filters ensure permanent pressure compensation and provide for optimal drainage in the system.
- ✓ Due to controlled irrigation of the bladder with the pump ball integrated in the feeding tube, post-operative tamponade of the bladder can be avoided.
- ✓ A doghouse on the bag ensures safe and hygienic positioning of the drain tap.



Capacity	Tube length	Article number	PZN
4.000 ml	120 cm	04/02/4000/TUR/GHC	06189004



#### Pump ball

The specific interaction of the pump ball and the two tube clamps allows to flush out blood clots in the bladder and tube system.



#### Profiled tube

The octagonal inner profile of the feeding tube assures excellent stability, but at the same time it is also flexible. A high lumen guarantees a safe drainage into the drip chamber in any situation.



#### Drip chamber

An anti-reflux valve divides the drip chamber and reliably prevents a backflow of urine and thereby obviates the risk against an ascending infection.



#### Drain tap in the doghouse

The special design of the piston in the drain tap ensures safe and easy operation and prevents leakages caused by improper use of the piston.

**Packaging Unit:** 10 pieces  
**Packaging:** sterile

Pictures may differ from the original

**CARE FLOW UB 4000**  
Closed urine drainage system



**drainage system, 4.000 ml**

- ✓ Can be fixed safely and correctly to many different types of bed frames by means of the twin holder with fixing cord.
- ✓ Bacteria- and fluid-tight venting filters ensure permanent pressure compensation and provide for optimal drainage in the system.
- ✓ Clear graduation markings allow for reliable volume control.
- ✓ A doghouse on the bag ensures safe and hygienic positioning of the drain tap.
- ✓ T-drain tap.
- ✓ Kink resistant tube.



Capacity	Tube length	Article number	PZN
4.000 ml	120 cm	04/02/4000/120/GHC	06718939



**Staggered connector**  
The staggered connector provides a reliable connection with the catheter. A needle-free sampling of urine is possible at any time, due to the needle-free sample port, which prevents needlestick injuries.



**Profiled tube**  
The octagonal inner profile of the feeding tube assures excellent stability, but at the same time it is also flexible. A high lumen guarantees a safe drainage into the drip chamber in any situation.



**Drip chamber**  
An anti-reflux valve divides the drip chamber and reliably prevents a backflow of urine and thereby obviates the risk against an ascending infection.



**Drain tap in the doghouse**  
The special design of the piston in the drain tap ensures safe and easy operation and prevents leakages caused by improper use of the piston.

**Packaging Unit:** 20 pieces  
**Packaging:** sterile

Pictures may differ from the original

**CARE FLOW UB 2000**  
Closed urine drainage system



**urine drainage system, 2.000 ml**

- ✓ Can be fixed safely and correctly to many different types of bed frames by means of the twin holder with fixing cord.
- ✓ Bacteria- and fluid-tight venting filters ensure permanent pressure compensation and provide for optimal drainage in the system.
- ✓ Clear graduation markings allow for reliable volume control.
- ✓ A doghouse on the bag ensures safe and hygienic positioning of the drain tap.
- ✓ Incorporated T-drain tap.
- ✓ Kink resistant tube.



Capacity	Tube length	Article number	PZN
2.000 ml	120 cm	04/02/2000/120/GHC	06718922
2.000 ml	90 cm	04/02/2000/90/GHC	09151597



**Staggered connector**  
The staggered connector provides a reliable connection with the catheter. A needle-free sampling of urine is possible at any time, due to the needle-free sample port, which prevents needlestick injuries.



**Profiled tube**  
The octagonal inner profile of the feeding tube assures excellent stability, but at the same time it is also flexible. A high lumen guarantees a safe drainage into the drip chamber in any situation.



**Drip chamber**  
An anti-reflux valve in the drip chamber provides for reliable prevention of urine reflux and highly efficient protection against ascending migration of bacteria.



**Drain tap in the doghouse**  
The special design of the piston in the drain tap ensures safe and easy operation and prevents leakages caused by improper use of the piston.

**Packaging Unit:** 20 pieces  
**Packaging:** sterile

Pictures may differ from the original

**CARE FLOW UB 2000S**  
Closed urine drainage system, soft tube



**urine drainage system, 2.000 ml**

- ✓ Urine drainage system with drip chamber prevents the risk of ascending germs and thus infections.
- ✓ Tube with large inner lumen made of soft material ensures very good drainage of the urine and high comfort of use.
- ✓ Twin holder with additional fastening cord enables secure attachment to different bed rails.
- ✓ Bacteria and liquid-tight ventilation filters ensure permanent pressure equalization.
- ✓ Clear scaling.
- ✓ Safe and at the same time hygienic positioning of the drain tap in the doghouse.
- ✓ Incorporated T-drain tap.



Capacity	Tube length	Article number	PZN
2.000 ml	120 cm	04/02/2000S/120	18728133
2.000 ml	90 cm	04/02/2000S/90	18728179

**Packaging unit:** 50 pieces  
**Packaging:** sterile

The **CARE FLOW BAG** serves as a sight protection for the **CARE FLOW UB 2000/UB 2000S** as discreet storage for the bed bags.

**CARE FLOW BAG**

- ✓ Handy Velcro fasteners make it quick and easy to attach the sight protection to the urine bag. A further Velcro fastener on the front allows a clean and uncomplicated usage of the drain tap.
- ✓ The sight protection is made of 100% polypropylene (PP) and is therefore resistant to organic acids, among other things.
- ✓ The sight protection is light gray, washable up to 40° and therefore suitable for reuse.



**Packaging Unit:** 10 pieces  
**Packaging:** non-sterile

Article number	PZN
04/99/01/945/GHC	09764331

Pictures may differ from the original

**CARE FLOW UB 2000A Disposable bedside bag**  
Sterile bedside bag with integrated twin holder



**bedside bag, 2.000 ml**

- ✓ Can be fixed safely and correctly to many different types of bed frames by means of the twin holder.
- ✓ An anti-reflux valve provides for reliable prevention of urine reflux and efficient protection against ascending migration of bacteria.
- ✓ Clear graduation markings allow for reliable volume control.
- ✓ A slit ensures the safe and hygienic positioning of the drain tap.
- ✓ T-drain tap.
- ✓ Kink resistant tube



Capacity	Tube length	Article number	PZN
2.000 ml	110 cm	04/04/2000/110/GHC	08871289



**Staggered connector**

The staggered connector provides a reliable connection with the catheter. A needle-free sampling of urine is possible at any time, due to the needle-free sample port, which prevents needlestick injuries.



**Profiled tube**

The octagonal inner profile of the feeding tube assures excellent stability, but at the same time it is also flexible. A high lumen guarantees a safe drainage into the drip chamber in any situation.



**Venting filter**

A bacteria and liquid-tight venting filters ensure permanent pressure equalisation. This ensures an optimal drainage of the system.



**Drain tap**

The special design of the piston in the drain tap ensures safe and easy operation and prevents leakages caused by improper use of the piston.

**Packaging Unit:** 50 pieces  
**Packaging:** sterile

Pictures may differ from the original

**CARE FLOW UB 2000E, UB 2000US**  
Sterile / non-sterile single-use bedside bag with holder



**bedside bag, 2.000 ml**

- ✓ The system is equipped with a universal staggered cone and can be connected to all common catheters.
- ✓ The product must be changed daily.
- ✓ An anti-reflux valve provides for reliable prevention of urine reflux and efficient protection against ascending migration of bacteria.
- ✓ Use the twin holder shown below to attach the bag to many different types of bed frames.
- ✓ The urine bag is emptied safely by means of the T-valve.



Article	Capacity	Sterility	Tube length	Article number	PZN
UB 2000E	2.000 ml	sterile	90 cm	04/04/2000/90/GHC	16047274
UB 2000US	2.000 ml	non-sterile	90 cm	04/06/2000/90/GHC	12418549
UB 2000US	2.000 ml	non-sterile	130 cm	04/06/2000/130/GHC	08845487

UB 2000E: **Packaging Unit:** 100 pieces, **Packaging:** sterile  
 UB 2000US: **Packaging Unit:** 100 pieces, **Packaging:** non-sterile **Sub-packaging:** each 10 pieces

**Holder for UB 2000E and UB 2000US**

- ✓ Plastic hanger for bedside bags with twin holder. Equipped with two pins taking up single-use bedside bags hanger holes.

Article number	PZN
04/99/BHL/PVC/GHC	06718951

**Packaging Unit:** 100 pieces **Sub-packaging:** 10 pieces each  
**Packaging:** non-sterile



Pictures may differ from the original

**CARE FLOW UB 2000USo, UB 1500USo, SB 1500USo**  
Non-sterile single-use bedside bag/secretion bag without drain tap



**single-use, 1.500 und 2.000 ml**

- ✓ The system is equipped with a universal staggered cone and can be connected to all common catheters.
- ✓ The product must be changed daily.
- ✓ An anti-reflux valve provides for reliable prevention of urine reflux and efficient protection against ascending migration of bacteria. (except the SB 1500USo).
- ✓ Use the twin holder shown below to attach the bag to many different types of bed frames.



Article	Capacity	Tube length	Anti-reflux valve	Article number	PZN
UB 2000USo	2.000 ml	90 cm	yes	04/07/2000/90/GHC	08845501
UB 1500USo	1.500 ml	90 cm	yes	04/07/1500/90/01	16839791
SB 1500USo	1.500 ml	90 cm	no	04/07/1500/90/00	16839785

UB 2000USo: **Packaging Unit:** 150 pieces, **Packaging:** non-sterile **Sub-packaging:** each 10 pieces  
 UB 1500USo: **Packaging Unit:** 200 pieces, **Packaging:** non-sterile **Sub-packaging:** each 10 pieces  
 SB 1500USo: **Packaging Unit:** 200 pieces, **Packaging:** non-sterile **Sub-packaging:** each 10 pieces

**Holder for UB 2000USo, UB 1500USo and SB 1500USo**

- ✓ Plastic hanger for bedside bags with twin holder. Equipped with two pins taking up single-use bedside bags hanger holes.

Article number	PZN
04/99/BHL/PVC/GHC	06718951

**Packaging Unit:** 100 pieces **Sub-packaging:** 10 pieces each  
**Packaging:** non-sterile



Pictures may differ from the original



**drainage systems, 1.000 ml**

- ✓ Closed urine drainage system.
- ✓ Bacteria- and fluid-tight venting filters provide for permanent pressure compensation and optimal drainage in the system.
- ✓ A doghouse on the bag ensures safe and hygienic positioning of the drain tap.
- ✓ T-drain tap.
- ✓ Needle free sampling port.
- ✓ Clear gradvation marks.



Capacity	Tube length	Article number	PZN
1.000 ml	30 cm	04/03/1000/30/GHC	06718885
1.000 ml	60 cm	04/03/1000/60/GHC	06718891
1.000 ml	90 cm	04/03/1000/90/GHC	06718916



**Stepped connector**  
 The special design of the connector steps provides for a reliable connection with the ring profile in the funnel end of the catheter. Needle-free urine sampling avoids needle injuries.



**Fleece coating**  
 The backside of the leg bag is coated with sweat-absorbing fleece. This effectively prevents irritation of the skin and considerably improves user comfort.



**Drip chamber**  
 An anti-reflux valve in the drip chamber provides for reliable prevention of urine reflux and highly efficient protection against ascending migration of bacteria.



**Drain tap**  
 The special design of the piston in the drain tap ensures safe and easy operation and prevents leakages caused by improper use of the piston.

**Packaging Unit:** 30 pieces  
**Packaging:** sterile

**Fixation for leg bags on page 54-58.**

Pictures may differ from the original



**1-chamber leg bag, 750 ml**

- ✓ The leg bag with 10 cm and 30 cm tube can conveniently be worn on the thigh.
- ✓ The leg bag with 50 cm tube comes with an additional stepped connector in the package. So the tube can be cut to the individual length and the leg bag can conveniently be worn on the thigh or lower leg.
- ✓ An anti-reflux valve provides for reliable prevention of urine reflux and efficient protection against ascending migration of bacteria.

Capacity	Tube length	Article number	PZN	Sterility
750 ml	10 cm	04/01/750/10/GHC	06718862	sterile
750 ml	30 cm	04/01/750/30/GHC	09772046	sterile
750 ml	50 cm	04/01/750/50/GHC	06718879	sterile
750 ml	70 cm	04/01/750A/70/GHC	19438785	sterile

**Packaging Unit:** 10 cm: 75 pieces, 30 cm: 60 pieces, 50 cm: 50 pieces, 70 cm: 50 pieces

Capacity	Tube length	Article number	PZN	Sterility
750 ml	10 cm	04/01/750US/10/GHC	01972716	non-sterile
750 ml	30 cm	04/01/750US/30/GHC	12341427	non-sterile
750 ml	50 cm	04/01/750US/50/GHC	00226543	non-sterile

**Packaging Unit:** 10 cm: 120 pieces, 30 cm: 100 pieces, 50 cm: 100 pieces, Sub-Packaging: 10 pieces each



**Stepped connector**  
 The special design of the connector steps provides for a reliable connection with the ring profile in the funnel end of the catheter. Needle-free urine sampling avoids needle injuries.



**Doghouse**  
 A doghouse on the bag ensures safe and hygienic positioning of the drain tap.



**Fleece coating**  
 The backside of the leg bag is coated with sweat-absorbing fleece. This effectively prevents irritation of the skin and considerably improves user comfort.



**Drain tap**  
 The special design of the piston in the drain tap ensures safe and easy operation and prevents leakages caused by improper use of the piston.

**Packaging:** sterile / non-sterile

**Fixation for leg bags on page 54-58.**

Pictures may differ from the original

## CARE FLOW UB 750F + UB 750UF sterile / non-sterile 1-way chamber leg bag systems with flip valve for mobile users



### 1-chamber leg bag, 750 ml

- ✓ Especially users with limited hand dexterity or reduced power can safely and hygienically empty this leg bag, due to our newly developed flip valve.
- ✓ The leg bag with 50 cm tube comes with an additional stepped connector in the package. So the tube can be cut to the individual length and the leg bag can conveniently be worn on the thigh or lower leg.
- ✓ An anti-reflux valve provides for reliable prevention of urine reflux and efficient protection against ascending migration of bacteria.



Capacity	Tube length	Article number	PZN	Sterility
750 ml	50 cm	04/01/750F/50/GHC	12371500	sterile

**Packaging Unit:** 50 pieces

Capacity	Tube length	Article number	PZN	Sterility
750 ml	50 cm	04/01/750UF/50/GHC	12371517	non-sterile

**Packaging Unit:** 100 pieces  
Sub-Packaging: 10 pieces each



#### Stepped connector

The special design of the connector steps provides for a reliable connection with the ring profile in the funnel end of the catheter. Needle-free urine sampling avoids needle injuries.



#### Doghouse

A doghouse on the bag ensures safe and hygienic positioning of the drain tap.



#### Fleece coating

The backside of the leg bag is coated with sweat-absorbing fleece. This effectively prevents irritation of the skin and considerably improves user comfort.



#### Drainage tap

The ergonomically formed flat **flip valve** of the drainage tube can be used easily and intuitively. To drain the leg bag the lever of the valve is flipped downwards by 180°.

**Packaging:** sterile / non-sterile

**Fixation for leg bags on page 54-58.**

Pictures may differ from the original

## CARE FLOW UB 500 + UB 500US sterile / non-sterile 3-chamber leg bag system for mobile users



### 3-chamber leg bag, 500 ml

- ✓ The leg bag with 10 cm and 30 cm tube can conveniently be worn on the thigh.
- ✓ The leg bag with 50 cm tube comes with an additional stepped connector in the package. So the tube can be cut to the individual length and the leg bag can conveniently be worn on the thigh or lower leg.
- ✓ An anti-reflux valve provides for reliable prevention of urine reflux and efficient protection against ascending migration of bacteria.



Capacity	Tube length	Article number	PZN	Sterility
500 ml	10 cm	04/01/500/10/GHC	09084990	sterile
500 ml	30 cm	04/01/500/30/GHC	12341203	sterile
500 ml	50 cm	04/01/500/50/GHC	09085009	sterile

**Packaging unit:** 10 cm: 75 pieces, 30 cm: 60 pieces, 50 cm: 50 pieces

Capacity	Tube length	Article number	PZN	Sterility
500 ml	10 cm	04/01/500US/10/GHC	09085015	non-sterile
500 ml	30 cm	04/01/500US/30/GHC	12341410	non-sterile
500 ml	50 cm	04/01/500US/50/GHC	09085021	non-sterile

**Packaging unit:** 10 cm: 120 pieces, 30 cm: 100 pieces, 50 cm: 100 pieces, Sub-Packaging: 10 pieces each



#### Stepped connector

The special design of the connector steps provides for a reliable connection with the ring profile in the funnel end of the catheter. Needle-free urine sampling avoids needle injuries.



#### 3-chambers

Due to the 3-chambers will the urine spread more regular through the system. In that way is it possible to wear the filled bag even more discrete beneath the clothing.



#### Fleece coating

The backside of the leg bag is coated with sweat-absorbing fleece. This effectively prevents irritation of the skin and considerably improves user comfort.



#### Drainage tap

The special design of the piston in the drain tap ensures safe and easy operation and prevents leakages caused by improper use of the piston.

**Packaging:** sterile / non-sterile

**Fixation for leg bags on page 54-58.**

Pictures may differ from the original

## CARE FLOW KB 1000 Closed urine drainage system for children



### Drainage System, 1.000 ml

- ✓ Closed urine drainage system.
- ✓ The bacteria- and liquid-tight ventilation filters ensure continuous pressure equalization, achieving optimal drainage.
- ✓ Twin holders and an additional cord allow for secure attachment to bed rails.
- ✓ A return tab guarantees safe and hygienic positioning of the drainage valve.
- ✓ Kink-resistant tubing.
- ✓ T-drain valve.
- ✓ Urine collection point.
- ✓ Clear volume scale.



Capacity	Tube length	Article Number	PZN
1.000 ml	90 cm	04/02/1000/90/GHC	19146953



**Staggered connector**  
The special alignment of the connector steps ensures a reliable connection with the reinforced ring at the catheter funnel. The needle-free urine sampling prevents needle-stick injuries.



**Twin holders**  
Twin holders are securely welded to the bag material to ensure a firm attachment to bed rails. Additional fixation can be achieved with the cord.



**Drip chamber**  
A reflux valve divides the drip chamber and reliably prevents urine backflow, effectively preventing the risk of ascending infections.



**Drain tap**  
A special counter-bulge on the piston of the drain valve guarantees easy and safe handling. Leaks caused by over-advancing the piston are prevented.

**Packaging unit:** 20 pieces  
**Packaging:** sterile

Pictures may differ from the original

## CARE FLOW Single-Chamber Urine Bag and CARE FIX Leg Bag Velcro Strap Set for Children



### Urine Bag, 210 ml

- ✓ Small urine bag for children and infants.
- ✓ The stepped connector ensures a reliable connection to the catheter and proper urine collection in the bag, even when positioned on lying children or infants.
- ✓ A reflux valve effectively prevents urine backflow from the bag.
- ✓ T-drain valve.
- ✓ Clear volume scale.



Volume	Article Number	PZN
210 ml	04/01/210/GHC	19146976

**Packaging unit:** 100 pieces  
**Packaging:** sterile

### Velcro Strap Set, 57 cm

- ✓ Comfortable-to-wear Velcro strap set for the UB 210 urine bag.
- ✓ The elastic material minimizes any disruption to blood circulation.
- ✓ 57 cm long Velcro straps, individually adjustable.

- 60°
- Do not tumble dry
- Do not iron

Length	Article Number	PZN
57 cm	05/01/BKL/KIN/GHC	19624393

**Packaging unit:** 10 sets, each with 2 straps  
**Packaging:** Non-sterile!



Pictures may differ from the original



### Universal Velcro

- ✓ Universal leg bag holding band with two attached Velcro straps to go through the upper eyelets of the leg bag and a special sewn on fixation for the tube of the leg bag system.
- ✓ The elastic material ensures that disturbances of the blood circulation in the thigh and lower leg are prevented as far as possible.
- ✓ The Velcro strap is 70 cm long and can be cut to the individual length. So it can be used on the thigh and lower leg.
- ✓ This universal Velcro strap can be used for the following urine leg bag systems: UB 1000, UB 750, UB 750F, UB 750US, UB 750UF, UB 500 und UB 500US.



- 60°
- Do not put in the dryer
- Do not iron

Lenght	Article number	PZN
70 cm	05/01/BKL/Uni/GHC	06718968

**Packaging Unit:** 10 pieces  
**Packaging:** sterile

Pictures may differ from the original



### Velcro strap premium

- ✓ sets consist of:  
an extra wide Velcro strap with two sewn-on Velcro loops to hold the punched eyelets on the upper part of the leg bag
- ✓ Both sets are equipped sewn-on loop to fixate the leg bag tube.
- ✓ The second Velcro Strap ensures a safe fixation of the lower part of the leg bag.
- ✓ Silicone strips on the back for a particularly good hold on the leg.
- ✓ The elastic material prevents a disruption of the blood circulation on the thigh or lower leg.
- ✓ 70 cm long Velcro tape can be shortened individually.
- ✓ The attached loop for an easier opening of the holding band.
- ✓ The wide premium velcro strap is also available separately, without the narrow hook and loop tape



#### Velcro strap with Velcro loops

Lenght	Article number	PZN
70 cm	05/01/BKL/PRE/SET	09241099
70 cm	05/01/BKL/PRE/GHC	06906196

**Packaging Unit:** 10 pieces  
**Packaging:** non-sterile

- 60°
- Do not put in the dryer
- Do not iron

- ✓ This Velcro strap premium set can be used for the following urine leg bag systems: UB 1000, UB 750, UB 750F, UB 750US, UB 750UF, UB 500 und UB 500US.
- ✓ The premium set with buttons (Art.-Nr. 05/01/BKL/PRE/03) comes with additional four buttons to fixate the upper part of the leg bag.

#### Velcro strap with buttons

Lenght	Article number	PZN
70 cm	05/01/BKL/PRE/03	12357931

**Packaging Unit:** 10 pieces  
**Packaging:** non-sterile



Pictures may differ from the original



**fixation stockings leg**

- ✓ Soft and ribbed lag pants with integrated leg bag pocket.
- ✓ Comfortable and safe fixation on the High as well as the lower leg.
- ✓ Due to the special design of the drain tap on the leg bag, emptying the bag is very simple.
- ✓ The extra-wide lateral Velcro fastener makes it very easy for the user to put on the leg pan (bei Art.-Nr. 05/02/10).
- ✓ The CARE FIX leg pants can be used for the following urine leg bag systems:  
 UB 750, UB 750F, UB 750US, UB 750UF, UB 500 und UB 500US.

**Fixation stockings with Velcro fastener**

Size	Article number	PZN
XXL	05/02/10/OS/XXL/GHC	08880957
XL	05/02/10/OS/XL/GHC	08871272
L	05/02/10/OS/L/GHC	08876513
M	05/02/10/OS/M/GHC	08845748
S	05/02/10/OS/S/GHC	07126649

**Measurement:** XXL = 60-85 cm, XL = 50-60 cm, L = 40-60 cm, M = 35-45 cm, S = 25-40 cm  
**Packaging Unit:** 10 pieces | **Packaging:** non-sterile

- 60° 60°
- Do not put in the dryer
- Do not iron

Measure in the middle of the upper or lower leg.

**Fixation stockings**

Size	Article number	PZN
XL	05/02/11/OS/XL/GHC	12418555
L	05/02/11/OS/L/GHC	07126017
M	05/02/11/OS/M/GHC	07126023
S	05/02/11/US/S/GHC	07126313

**Measurement:** XL = 65-95 cm, L = 55-80 cm, M = 45-65 cm, S = 30-45 cm  
**Packaging Unit:** 10 pieces | **Packaging:** non-sterile



**Upper and lower leg cuffs**

- ✓ The leg cuffs are made of soft, washable cotton blend fabric with an integrated leg bag pocket.
- ✓ Silicone strips on the back for a particularly good hold on the leg.
- ✓ Wearable on both sides.
- ✓ Leg bags are held in the integrated pocket that can be opened with Velcro, including opening for drain tap.
- ✓ The CARE FIX cuffs can be used for the following urine leg bag systems:  
 UB 1000, UB 750, UB 750A, UB 750F, UB 750US, UB 750UF, UB 500 und UB 500US.

It is measured in the middle of the thigh.

**Leg bag holder for the thigh**

Size	Article number	PZN
XL	05/02/12/XL/GHC	09080408
L	05/02/12/L/GHC	09080383
M	05/02/12/M/GHC	09080377

**Measurement:** XL = 61 cm, L = 58 cm, M = 54 cm  
**Packaging Unit:** 10 pieces  
**Packaging:** non-sterile

- 60° 60°
- Do not put in the dryer
- Do not iron

It is measured in the middle of the lower leg.

**Lower Leg Cuffs**

Size	Article number	PZN
L	05/02/13/L/GHC	09080437
M	05/02/13/M/GHC	09080420
S	05/02/13/S/GHC	09080414

**Measurement:** L = 37 cm, M = 35 cm, S = 32 cm  
**Packaging Unit:** 10 pieces  
**Packaging:** non-sterile



**CARE FIX**  
Premium one leg pants



**Premium one leg pants**

- ✓ The premium one leg pants are made from a washable cotton blend has an integrated leg bag pocket.
- ✓ Individually adjustable belt with silicone strips for a safe fixation.
- ✓ The premium one leg pant can be worn on both the right as well as the left thigh.
- ✓ Leg bags are positioned into the integrated pocket with a Velcro fastener; including an opening for the drain tap.
- ✓ Easy emptying of the leg bags due to the comfortable positioning of the drain tap.
- ✓ The leg pants can be used for the following urine leg bag systems: UB 1000, UB 750, UB 750A, UB 750F, UB 750US, UB 750UF, UB 500 und UB 500US.



- 60°
- Do not put in the dryer
- Do not iron

Measure in the middle of the upper or lower leg.

Size	Article number	PZN
XL	05/03/10/XL/GHC	09080472
L	05/03/10/L/GHC	09080466
M	05/03/10/M/GHC	09080443

**Measurement:** XL = 98-120cm, L = 90-100cm, M = 85-95cm

**Packaging Unit:** 10 pieces

**Packaging:** non-sterile

Pictures may differ from the original

**CARE DRESS**  
Suprapubic dressing set



**Dressing set**

- ✓ Suprapubic dressing set for optimal care of the suprapubic exit site.
- ✓ The included smaller plaster serves to temporarily fix the suprapubic catheter to the abdominal wall, thus preventing the catheter from shifting.



**Included in this set:**

- 1 x elastic plaster 10 x 10 cm
- 1 x elastic plaster 5 x 3 cm
- 3 x non-woven swabs 5 x 5 cm, 6-ply
- 2 x non-woven split swabs 5 x 5 cm, 6-ply
- Packed in a peel-apart pouch -

Article number	PZN
08/02/VS/GHC	06718974

**Packaging Unit:** 25 sets

**Packaging:** sterile

Pictures may differ from the original



### Dressing set

- ✓ Abdominal PEG dressing set for optimal care of the PEG exit site.



#### Included in this set:

- 1 x elastic plaster 12 x 15 cm
- 2 x non-woven swabs 7,5 x 7,5 cm, 4 ply
- 2 x non-woven split swabs 7,5 x 7,5 cm, 4 ply
- Packed in a peel-apart pouch -

Article number PZN  
08/02/PEG/GHC 05461059

**Packaging Unit:** 15 sets  
**Packaging:** sterile



CARECATH *Liv&Leiv*

Optimum length for safe bladder emptying.



- ✓ Excellent lubricity.
- ✓ Atraumatic, rounded drainage eyelets.
- ✓ Ready-to-use.
- ✓ Touch-free application.
- ✓ Easy to handle and guide.

Pictures may differ from the original



**IC disposable catheter**

- ✓ Very good gliding properties thanks to hydrophilic coating.
- ✓ Ready for immediate use thanks to ready-to-use, hydrophilic coating.
- ✓ Atraumatic application due to rounded drainage eyes.
- ✓ Optimal length of 17 cm ensures complete bladder emptying.
- ✓ Non-contact application and good guideability due to the correct rigidity of the catheter.
- ✓ Good grip of the handle allows good handling, even with limited dexterity.
- ✓ Easy connection of urine collection bag and/or instillation solution via the handle cone.



Atraumatic, rounded drainage eyes



Grip of the handle

CARE CATH *Leiv*

Size	Length	Article number	PZN
CH 10	17 cm	LIV-H1017	18890885
CH 12	17 cm	LIV-H1217	18890916
CH 14	17 cm	LIV-H1417	18890922
CH 16	17 cm	LIV-H1617	18890939

**Packaging unit:** 30 pieces  
**Packaging:** sterile

Pictures may differ from the original



**IC disposable catheter**

- ✓ Very good gliding properties thanks to hydrophilic coating.
- ✓ Ready for immediate use thanks to ready-to-use, hydrophilic coating.
- ✓ Atraumatic application due to rounded drainage eyes.
- ✓ Non-contact application and good guideability due to the correct rigidity of the catheter.
- ✓ Good grip of the handle allows good handling, even with limited dexterity.
- ✓ Easy connection of urine collection bag and/or instillation solution via the handle cone.



Atraumatic, rounded drainage eyes



Grip of the handle

CARE CATH *Leiv*

Size	Length	Article number	PZN
CH 12	41 cm	LEIV-H1241	18890833
CH 14	41 cm	LEIV-H1441	18890856
CH 16	41 cm	LEIV-H1641	18890862

**Packaging Unit:** 30 pieces  
**Packaging:** sterile

Pictures may differ from the original

## Curan Lady & Man IC disposable catheter with hydrogel coating, for woman & man



### IC disposable catheters

- ✓ Compact disposable catheter for discreet use.
- ✓ Good gliding properties thanks to its hydrogel coating.
- ✓ Polished drainage eyes enable comfortable insertion and prevent irritation or injury to the urethral mucosa.
- ✓ Disposable catheters are ready to use.
- ✓ The catheter material is phthalate-free.
- ✓ Funnel offers connection option for urine collection bag or bladder instillation.



**Curan Lady**  
Catheter length: 12 cm



**Curan Man**  
Catheter length: 30 cm

#### Curan Lady

Size	Length	Article number	PZN
CH 08	12 cm	CURAN-CL08	19341962
CH 10	12 cm	CURAN-CL10	19341979
CH 12	12 cm	CURAN-CL12	19341985
CH 14	12 cm	CURAN-CL14	19341991

**Packaging Unit:** 30 pieces  
**Packaging:** sterile

#### Curan Man

Size	Length	Article number	PZN
CH 12	30 cm	CURAN-CM12	19342022
CH 14	30 cm	CURAN-CM14	19342039
CH 16	30 cm	CURAN-CM16	19342045

**Packaging Unit:** 30 pieces  
**Packaging:** sterile



Atraumatic catheter eyes

## Curan Complete & Closed System IC disposable catheter with integrated collection bag, for women & men




### IC disposable catheter

- ✓ Complete systems with disposable catheter and integrated collection bag.
- ✓ Good gliding properties thanks to hydrogel coating.
- ✓ Polished drainage eyes enable comfortable insertion and prevent irritation or injury to the urethral mucosa.
- ✓ Disposable catheters are ready to use.
- ✓ The catheter material is phthalate-free.
- ✓ The introducer tip can be inserted hygienically into the urethra to bypass the germ-colonised urethral orifice by not coming into contact with the catheter itself.



**Curan Complete compact system**  
Collection bag: 1.000 ml

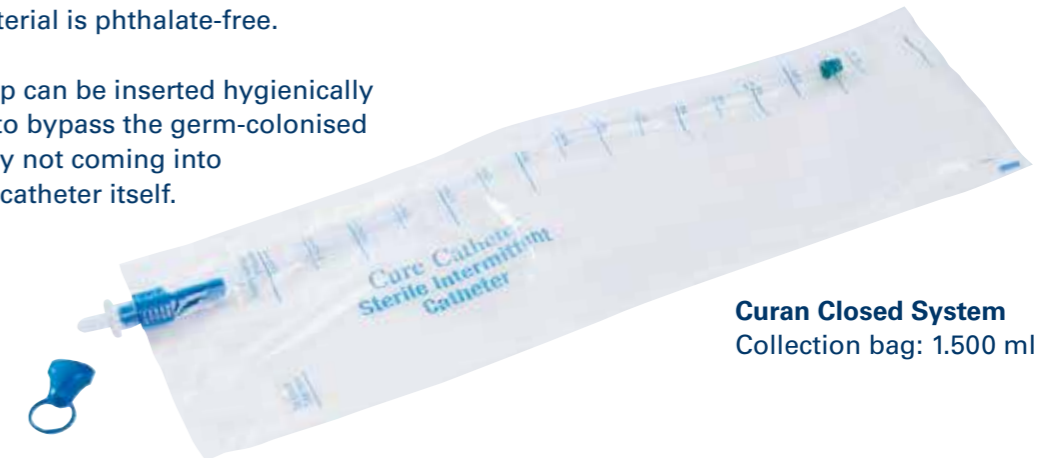

**Curan Closed System**  
Collection bag: 1.500 ml

#### Curan Complete with collection bag, 1.000ml

Size	Length	Article number	PZN
CH 12	40 cm	CURAN-CC12	19341933
CH 14	40 cm	CURAN-CC14	19341956

**Packaging Unit:** 30 pieces  
**Packaging:** sterile

#### Curan Closed System with collection bag, 1.500ml

Size	Length	Article number	PZN
CH 10	40 cm	CURAN-CB10	19341850
CH 12	40 cm	CURAN-CB12	19341867
CH 14	40 cm	CURAN-CB14	19341896
CH 16	40 cm	CURAN-CB16	19341904

**Packaging Unit:** 25 pieces  
**Packaging:** sterile



Introducer tip

Pictures may differ from the original

Pictures may differ from the original



### Safety port cannula

- ✓ Handling and insertion into the port is easy thanks to the very well and firmly manufactured gripper with integrated grooves.
- ✓ Thanks to the excellent Huber cut, the port cannula can be inserted into the port without coring and with virtually no pain. As a result, no skin or silicone particles are introduced into the bloodstream.
- ✓ The clamp can be firmly locked and released again. The colour-coded tube clamps make it easy to clearly distinguish between the sizes.
- ✓ When the port cannula is removed, the needle locks so firmly into the base plate that no needlestick injuries and infections can occur. The requirements of TRBA 250 are met and health-care professionals are protected.

Size	Lenght	Article number	PZN
18G	19 mm	09/01/18G/19	19223381
18G	25 mm	09/01/18G/25	19223398
19G	12,5 mm	09/01/19G/12	19707274
19G	15 mm	09/01/19G/15	19146203
19G	19 mm	09/01/19G/19	19146226
19G	25 mm	09/01/19G/25	19146255
19G	32 mm	09/01/19G/32	19146261
19G	38 mm	09/01/19G/38	19146278

18G available in the course of 2025.

Size	Lenght	Article number	PZN
20G	12,5 mm	09/01/20G/12	19707245
20G	15 mm	09/01/20G/15	19146284
20G	19 mm	09/01/20G/19	19146290
20G	25 mm	09/01/20G/25	19146309
20G	32 mm	09/01/20G/32	19146315
20G	38 mm	09/01/20G/38	19146321
22G	12,5 mm	09/01/22G/12	19707268
22G	15 mm	09/01/22G/15	19146338
22G	19 mm	09/01/22G/19	19146344
22G	25 mm	09/01/22G/25	19146350
22G	32 mm	09/01/22G/32	19146367
22G	38 mm	09/01/22G/38	19146373

**Packaging Unit:** 15 pieces  
**Packaging:** sterile



Good handling and Easy grip due to wingless gripper with integrated grooves.



Protection against needlestick injuries, due to safety mechanism. The needle locks into the base plate after the removal of the port cannula.



Coring-free application thanks to sharp Huber cut and good wearing comfort due to flat contact surface.



Firm closure of the connecting tube by means of a tube clamp, which supports therapy and safety.



Visit our website [www.ghc-medical.de](http://www.ghc-medical.de). There you will find detailed information about our company as well as the entire product portfolio.

Our sales team is looking forward to your call:

Mon.-Thurs.: 08:00 – 16:00 Uhr

Fri.: 08:00 – 14:00 Uhr

**GHC** German  
Health Care GmbH



**GHC German Health Care GmbH**

Birklück 15  
24999 Wees/Flensburg

Tel. +49 (0) 4631 / 572 00 0  
Fax +49 (0) 4631 / 572 00 19

info@ghc-medical.de  
www.ghc-medical.de

Version 06/2025