

Ubiq[®]

Therapeutic Comfort Shoe

NEW



Ubiq[®] Way



**The
everyday hero.
Style included.**



More lightness in every step.

The therapeutic comfort shoe Ubiq® is designed for people who value comfortable, dynamic walking. Its natural, smooth rolling motion ensures an all-around pleasant walking experience in everyday life. Functional, stylish, and so versatile that it fits anywhere – whether at work, during leisure time, or at home.

The name Ubiq® is derived from the word ubiquitous, meaning “present everywhere”. That’s exactly what this shoe represents: Universal acceptance, flexibility, and versatility. A true allrounder that fits every situation and complements every style.



dynamic
shock-absorbing
non-slip



Indications

- Mild to moderate foot / toe deformities such as hallux valgus, hammer toes or claw toes
- Hallux limitus
- Diabetes mellitus without PNP / PAOD (risk group 0 according to the German Diabetes Society)
- Diabetes mellitus without PNP / PAOD with foot deformity (risk group 1 according to the German Diabetes Association)

Features



Outsole

- Back-shifted ball rocker
→ Forefoot off-loading
- Integrated, large, abrasion-resistant rubber coating (4 mm)
→ Non-slip
- Shock-absorbing, force-absorbing sole construction
- Continuous and reinforced rigid shoe base
→ Forefoot off-loading
- Reduced footprint
→ Dynamic, fast rolling motion
- Lower heel-to-toe drop (7 mm)
→ Natural foot position



Insole

- Anatomically moulded footbed (6 mm)
→ Moderate heel guide
→ Slight lengthwise and crosswise arch support
→ PU sole with smooth microfibre cover for good cushioning
+ easy foot insertion
- Replaceable
- Sufficient space (6 mm) for customised insoles



Shoe Last Ubiq®

- Foot-friendly last shape
- Ideal ball width
- Tailored shape
- Toe box height 23.3 – 27.6 mm (depending on size)



Lacing

- Elasticated laces with quick-release fastener
→ Easy to put on and take off
- Non-elastic, flat laces (spare set enclosed)
→ Better support



Shaft

- Air-permeable
- Reinforced toe
- Anatomical ankle cutout
→ Ideal rearfoot fit
- Optimally adapted heel cup
→ Substantial and cushioned heel support
- Loop on the heel
→ Easy to put on
- Stable upper material
- Soft liner
→ Good padding
→ Protection against pressure and external influences
→ Comfortable to wear

Black




Grey



Navy



Lemon

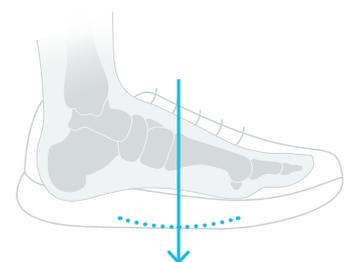


Back-shifted, dynamic ball rocker



Sizes + Measurements

Ubiq® Way

EU Shoe size	Insole		Foot length	Ball girth	Toe box height	Art. No.			
	A	B				Black	Grey	Navy	Lemon
36	24.0	9.0	22.5 – 23.1	22.5	2.33	UBWBA36	UBWGY36	UBWNY36	UBWLM36
37	24.6	9.1	23.1 – 23.7	23.0	2.37	UBWBA37	UBWGY37	UBWNY37	UBWLM37
38	25.3	9.2	23.8 – 24.4	23.4	2.41	UBWBA38	UBWGY38	UBWNY38	UBWLM38
39	26.0	9.3	24.5 – 25.1	23.8	2.45	UBWBA39	UBWGY39	UBWNY39	UBWLM39
40	26.6	9.4	25.1 – 25.7	24.3	2.49	UBWBA40	UBWGY40	UBWNY40	UBWLM40
41	27.3	9.6	25.8 – 26.4	24.7	2.52	UBWBA41	UBWGY41	UBWNY41	UBWLM41
42	28.0	9.8	26.5 – 27.1	25.2	2.56	UBWBA42	UBWGY42	UBWNY42	UBWLM42
43	28.6	10.1	27.1 – 27.7	25.6	2.60	UBWBA43	UBWGY43	UBWNY43	UBWLM43
44	29.3	10.2	27.8 – 28.4	26.1	2.64	UBWBA44	UBWGY44	UBWNY44	UBWLM44
45	30.0	10.3	28.5 – 29.1	26.5	2.68	UBWBA45	UBWGY45	UBWNY45	UBWLM45
46	30.6	10.4	29.1 – 29.7	27.0	2.72	UBWBA46	UBWGY46	UBWNY46	UBWLM46
47	31.3	10.5	29.8 – 30.4	27.4	2.76	UBWBA47	UBWGY47	UBWNY47	UBWLM47

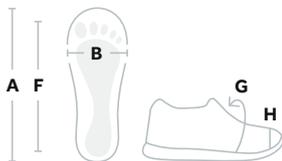

All measurements are in cm and refer to the insole length.
The indicated shoe sizes are for guidance only.
The actual size depends primarily on the volume of the bandage or dressing.



Find the right size for you



Exchangeable insole



For diabetics



Light weight



DARCO (Europe) GmbH

Gewerbegebiet 18 | 82399 Raisting | Germany | Tel. +49 8807 9228-0 | info@darco.de | www.darco.de

