





Intended purpose:	mediform is an orthopaedic medical device used to create a negative impression of feet.		
Function:	The professional user inserts the patient's foot into the FSG impression plate to create a precise impression of the foot. This precise impression keeps its shape, as the foam does not spring back after pressure has been applied to it.		
Sectors:	Insole production, custom shoe production, archiving, prostheses / orthoses, identification of prominent areas, determination of measurements		
Possible further processing:	Casting with plaster or resin, scanning, among other things.		
Product description:	There are two sheets of rigid impression foam encased in a collapsible box made of cardboard. These are then packed in the desired number into one of the shipping units available.		
Materials:	Foam = phenolic resin-based	Collapsible box = cardboard	Shipping box = cardboard
Specifications:	Variant Ideal 38		
Foam sheet dimensions:	330x145x38		
Foam sheet weight:	17g		
Collapsible box dimensions:	340x150x80		
Collapsible box weight:	102g		
Foam bulk density:	8-12 kg/m ³		
Available foam colours:	Pink, blue, green, yellow, violet, pink, red, grey		

Diagram with measurements:

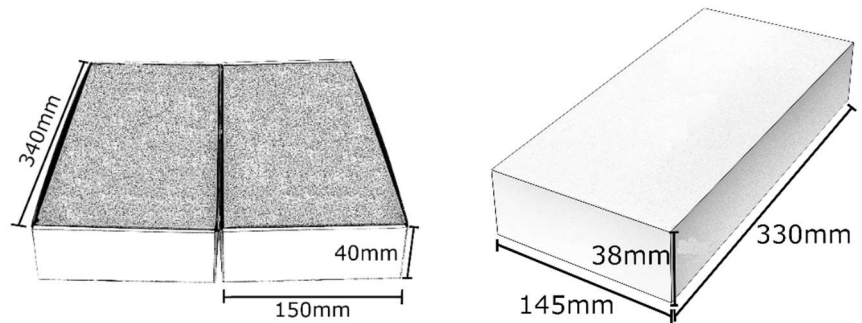

Available shipping units	25 pack	40 pack	50 pack
Weight	3.5 kg	5.5 kg	7 kg
Dimensions	77x42x35	70x60x42	77x70x42
Overall dimensions	231cm	274cm	301cm
Volume	0.12	0.17	0.23
			
Storage	Dry at temperatures between 18 and 24° C with a humidity between 55 and 65% and store shipping boxes upright. Do not store more than 4 boxes on top of each other.		
Disposal	FSG imprint foam can be disposed of in household waste. Waste codes 20 01 39 (municipal waste as well as commercial, made from plastic) The cardboard box can be added to waste paper recycling 15 01 01 (packaging made from paper and cardboard)		

FSG

Formschaum Gronau GmbH

Marie-Curie-Str.18

D-48599 Gronau

Tel. +49 (0) 2562 907 676 0

Fax: +49 (0) 2562 907 676 6

Email: info@formschaum-gmbh.de

More information at formschaum-gmbh.de

Hinweise:

Einmalprodukt

Vor Feuchtigkeit schützen

Nicht werfen, Bruchgefahr

